



Modernizing legacy Vaccine processes

Dr. Mats Lundgren
GE Healthcare Life Sciences

Imagination at work

A Modern Solution for Acellular Pertussis Vaccine

123RF®

Whole-cell (wP) - Acellular Pertussis (aP)

wP Vaccines

70 year old technology based on killed *B. pertussis* strains

High protection efficiency ~78%

Associated with side effects and safety concerns

The reactogenicity of wP vaccine was thought to be too high to permit routine use in older children, adolescents and adults.

aP Vaccines

Introduced in 1990's

aP contain ≥ 1 of the separately purified antigens: pertussis toxin (PT), filamentous hemagglutinin (FHA), pertactin (PRN), and fimbriae (FIM) type 2 and 3.

aP is now the dominant type in the industrialized world

aP containing vaccines with reduced concentrations of the antigen have been formulated for use in adolescents and adults




Traditional Process vs. Modern Solution

Challenges

- Time consuming
- Unable to purify separate antigens
- Low purity
- Low yield
- Difficult to scale up
- Salt disposal issue

Advantages

- Able to purify separate antigens
- High yield
- High purity
- Easy to scale up
- Time saving
- Environmentally friendly

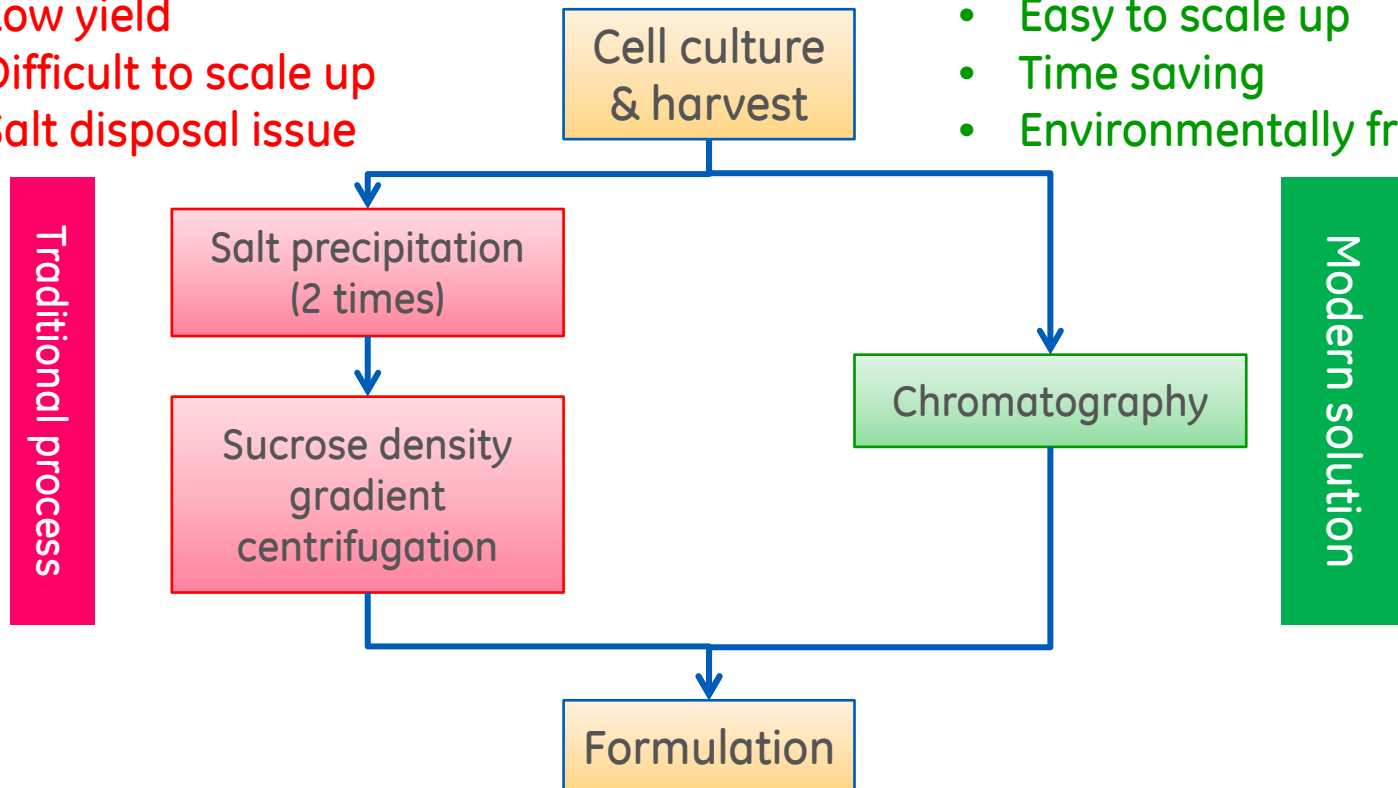

Project Goal

Traditional process

Chinese pharmacopeia requirement and current situation

- Contain 2 antigens:
Pertussis toxoid (PT),
Filamentous Hemagglutinin (FHA)
- Purity >85% (SDS-PAGE)

Yield around 10%

Lack of stable antigen quantitative assay

Current Project

Develop a modern process for pertussis vaccine


- Contain 3 antigens:
Pertussis toxoid (PT),
Filamentous Hemagglutinin (FHA)
Pertactin (PRN)
- Purity >95% (SDS-PAGE)

Yield >30%


Establish quantitative antigen determination using Biacore™ platform



Design overview

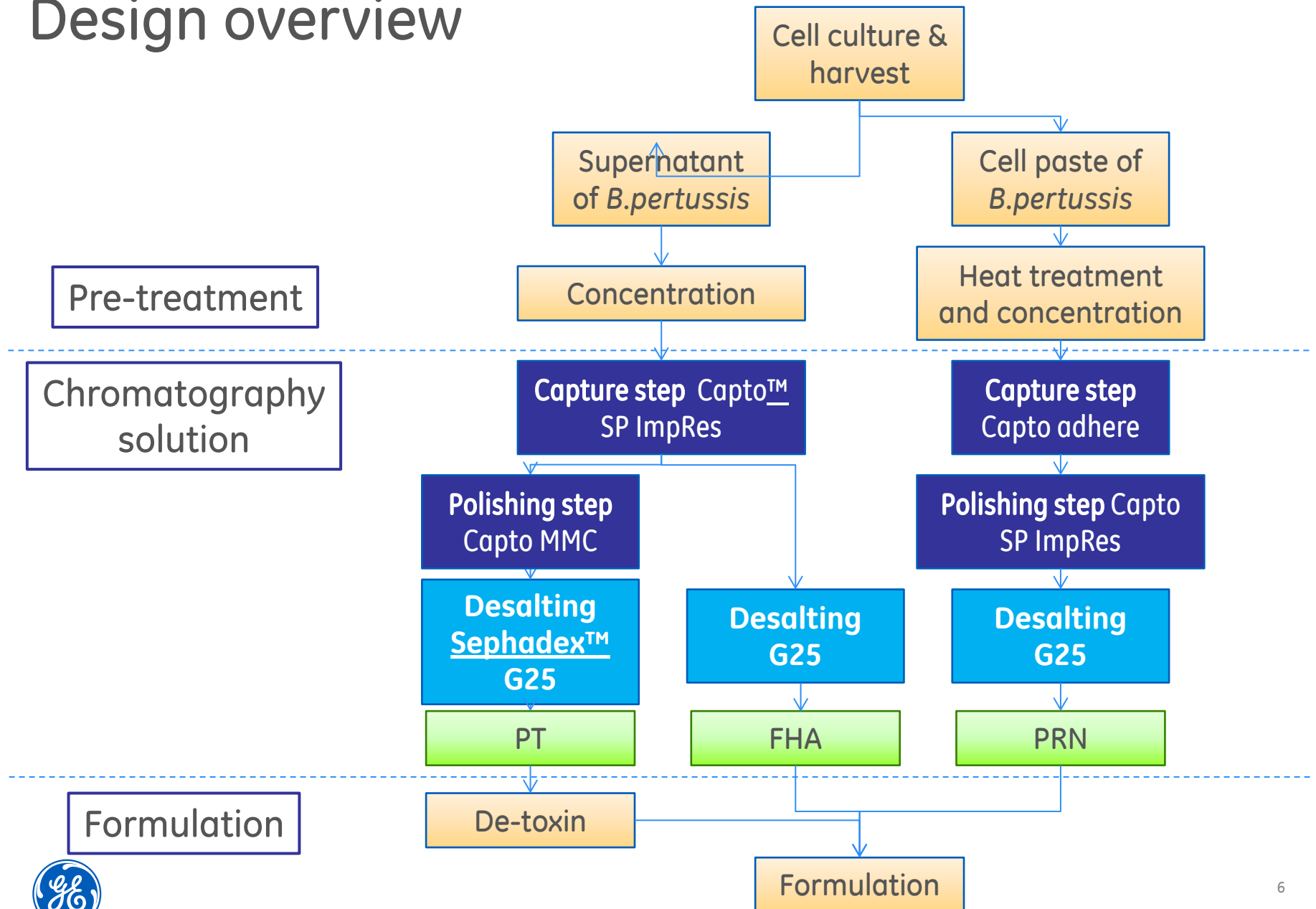

Design overview

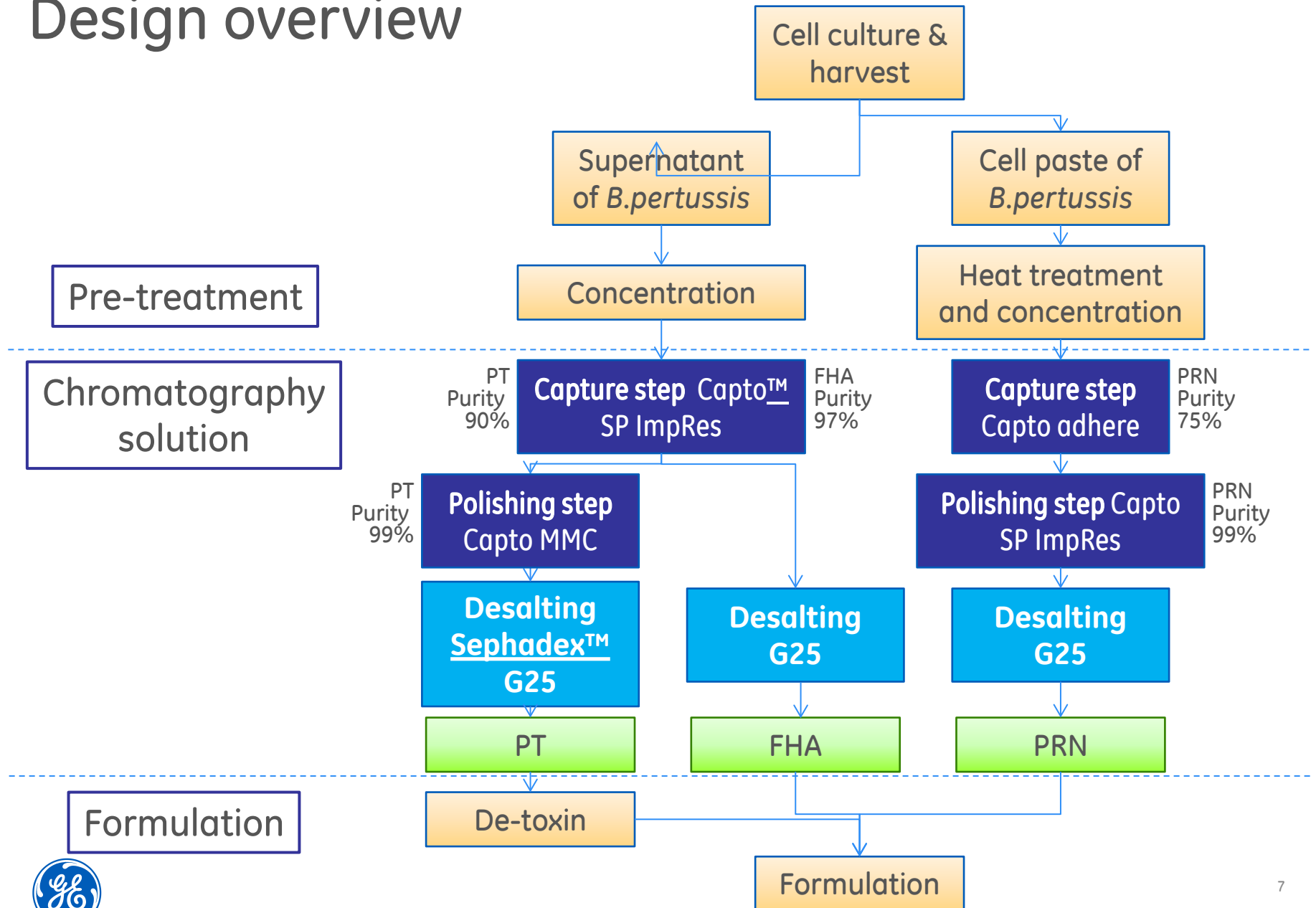

Process Highlights

1. Modern process to produce PT, FHA & PRN using bioprocess friendly, easily scalable, new generation chromatography platform.
2. Environmentally friendly.
3. Increase purity from 85% to >95%.
4. Reduce manufacture time from month to days.
5. Recovery increased from 10% to 30%.
6. Establish a sensitive, stable platform using Biacore to quantify PT & FHA.



Modern Process for Meningococcal Vaccine



Meningococcal Vaccine



13 clinically significant serotypes. A, B, C, W-135, Y responsible for 90% of global cases


Vaccine for A, C, W, Y are produced using capsular polysaccharide (PS), conjugant technology to enhance immunogenicity




Traditional Process vs. Modern Solution

Traditional process challenge:

- Phenol use
- High-speed centrifuge

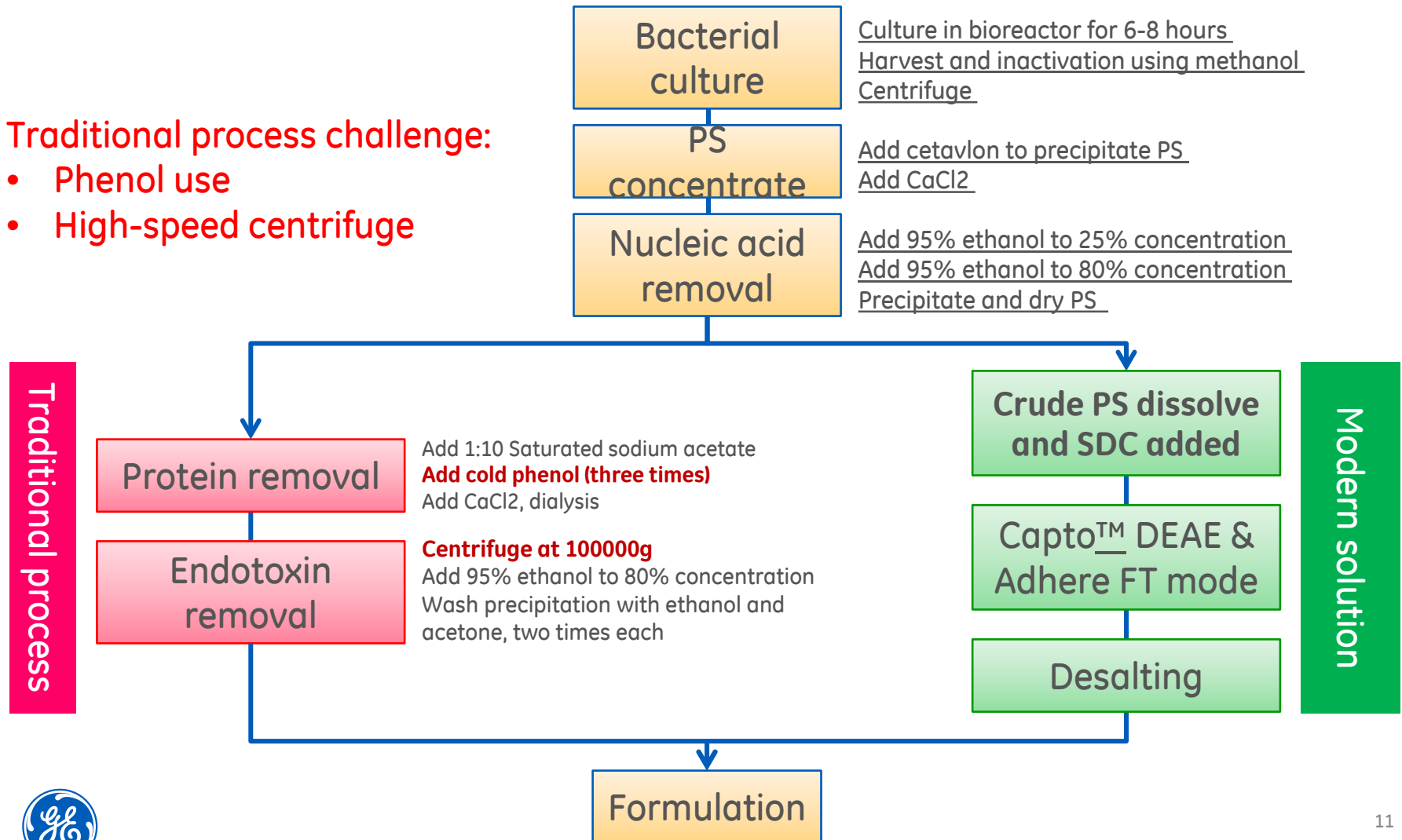


Typical Result From Upgraded Process 1x Capto DEAE + 2x Capto Adhere

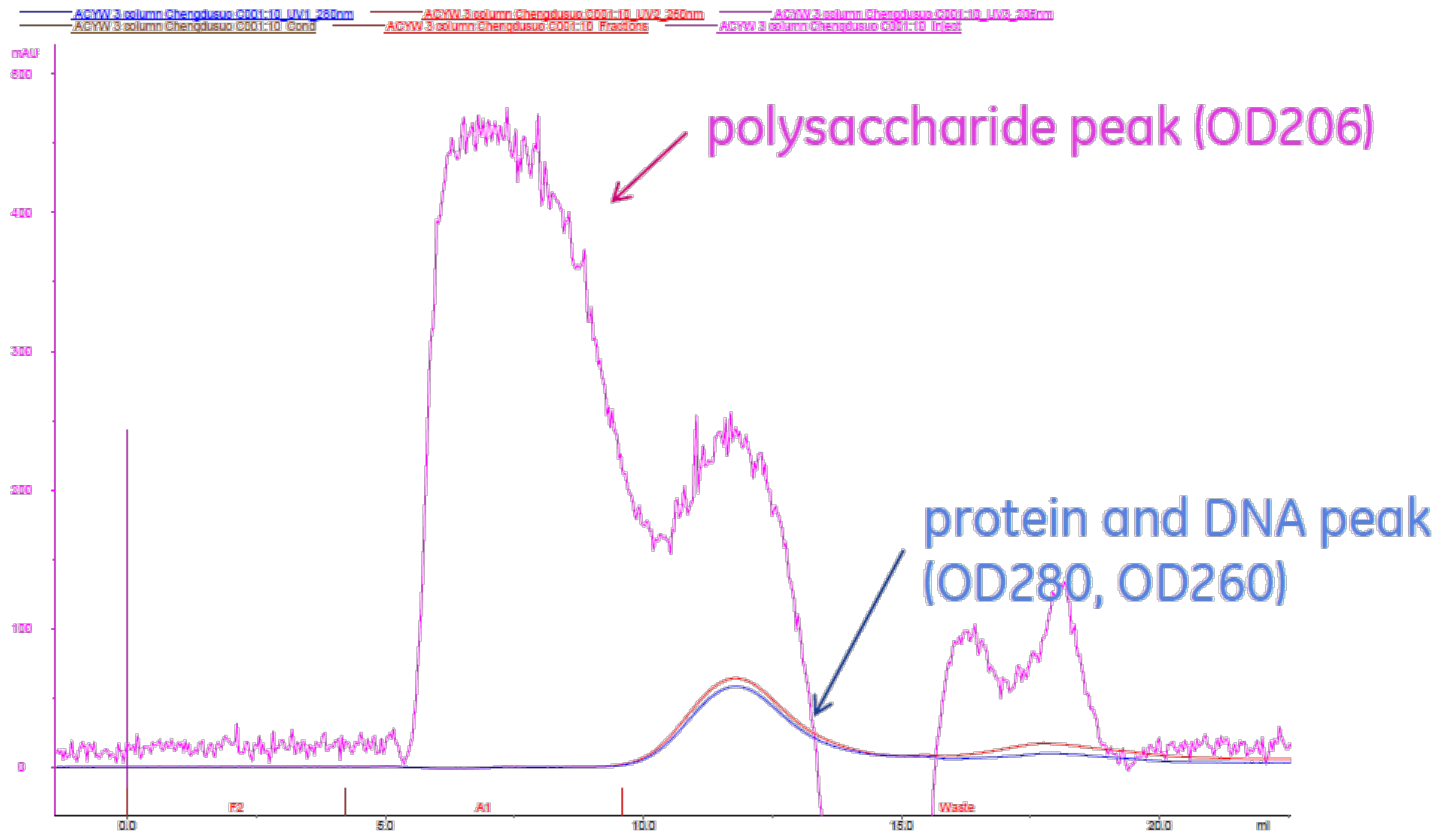

ALEX + HIC + SEC effect!




Separation Flowchart

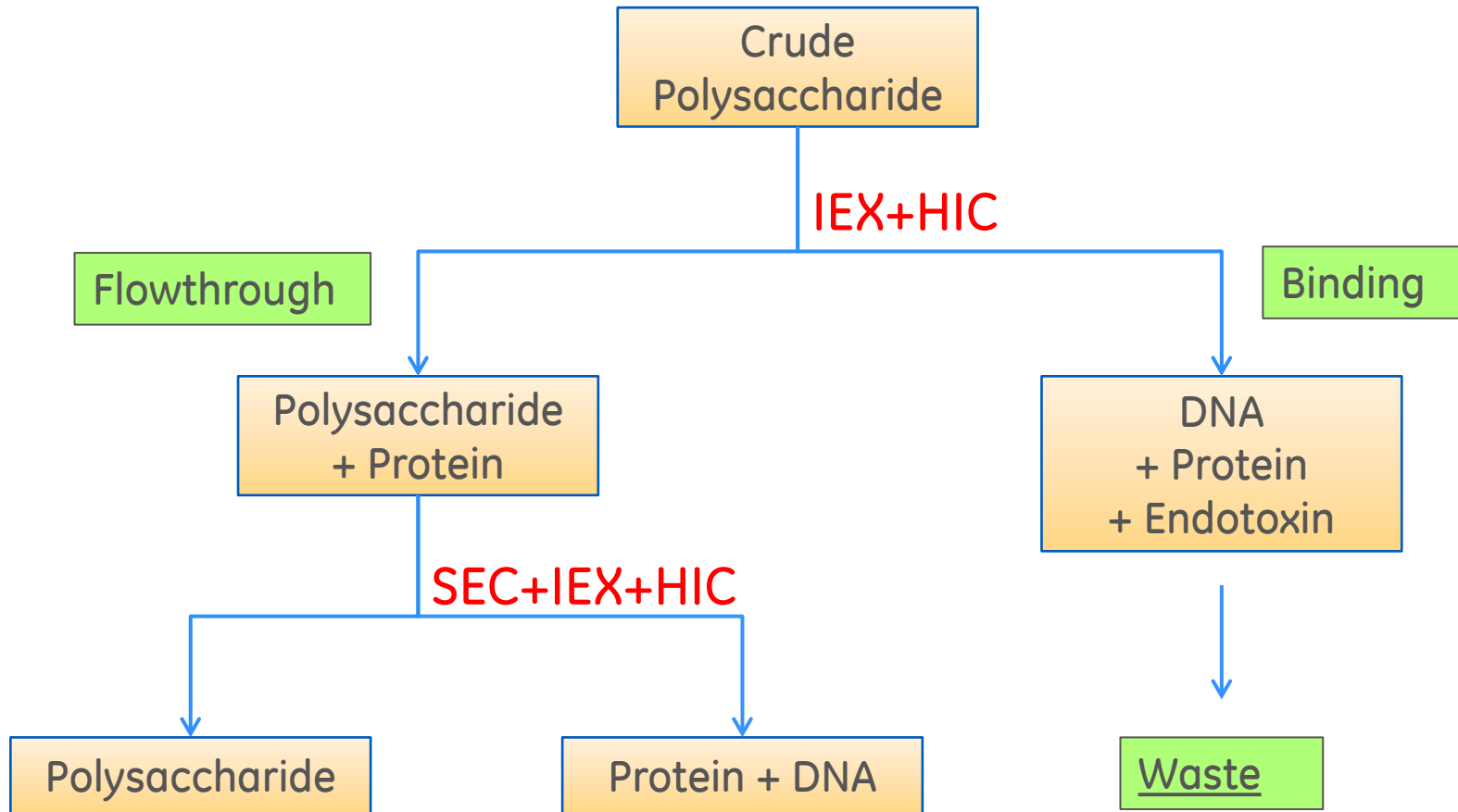

Modern solution for Meningococcal Vaccine A,C,W,Y

Advantages vs. traditional process:

- No phenol use in process, benefit environment & operator's health & safety
- Easy to scale up
- Simple flow-through mode
- All 4 serotypes using same process
- Protein/DNA/endotoxin in products meet requirement

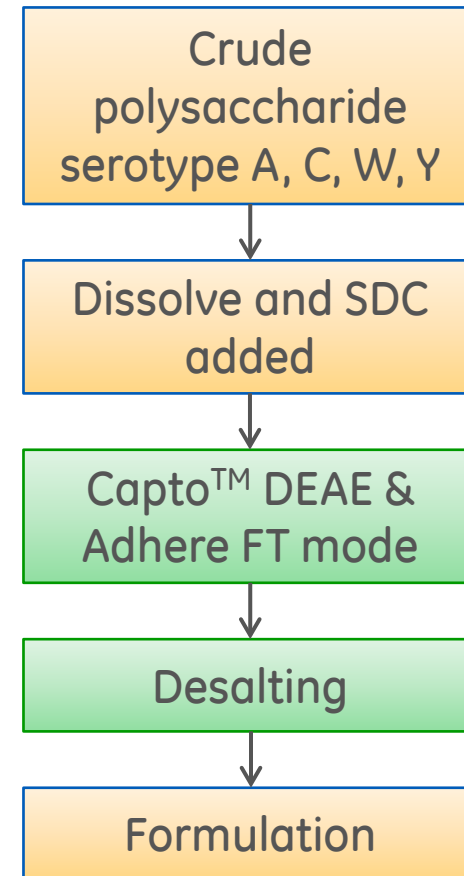

Live Influenza virus production



Virus seed stock safety

EP citation:

"Seed lots/cell banks. The master seed lot or cell bank is identified by historical records that include information on its origin and subsequent manipulation. Suitable measures are taken to ensure that no extraneous agent or undesirable substance is present in a master or working seed lot or a cell bank."



Influenza process overview



Scale-up from small scale to pilot scale in single-use format


Comparison of culture performance in 10 L and 50 L microcarrier culture in rocking bioreactors

Downstream purification in flow-through chromatography mode with Capto™ Q and Capto Core 700 chromatography media (resins)




Cell growth in single-use bioreactor stage

Cell concentration

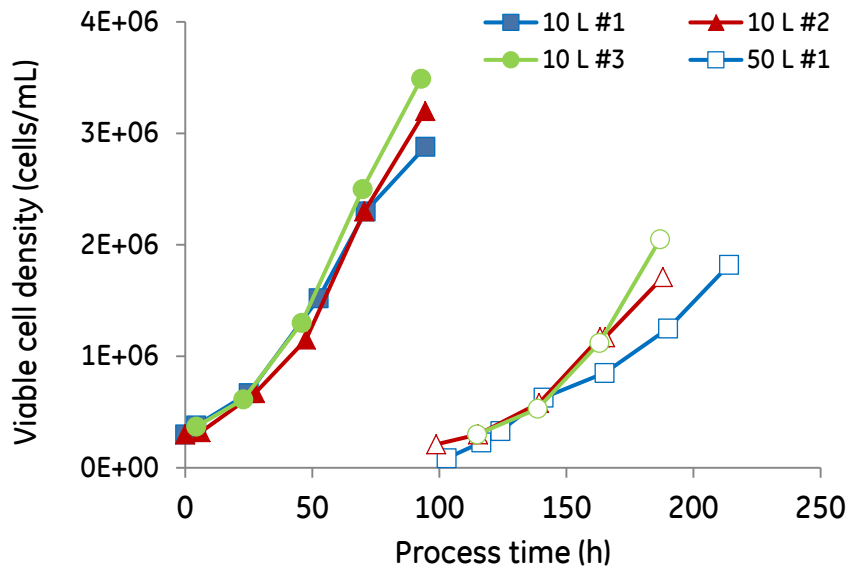

Average growth rate

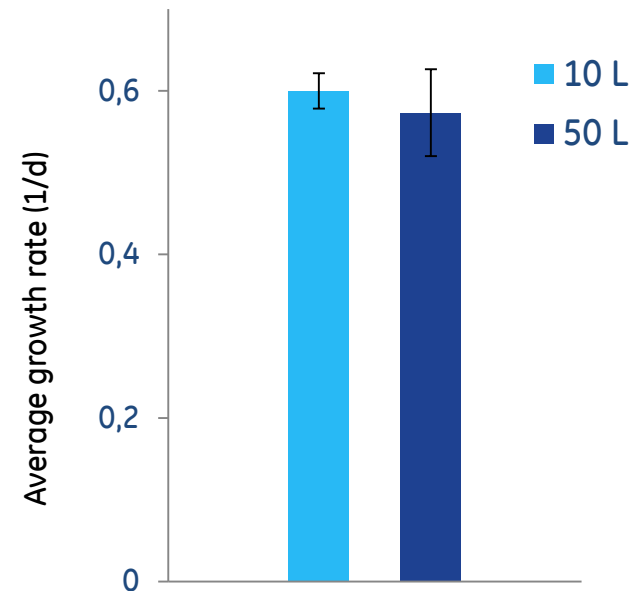


Bead to bead transfer was successful and cell growth was comparable at 10 L and 50 L scale




Virus growth kinetics

HA concentration and virus titer during culture

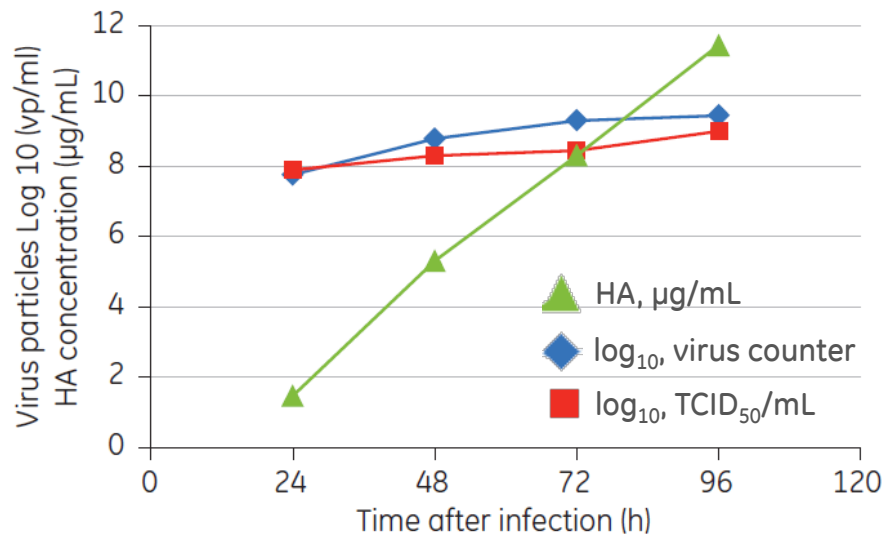


Cell morphology at time of harvest (96 h)

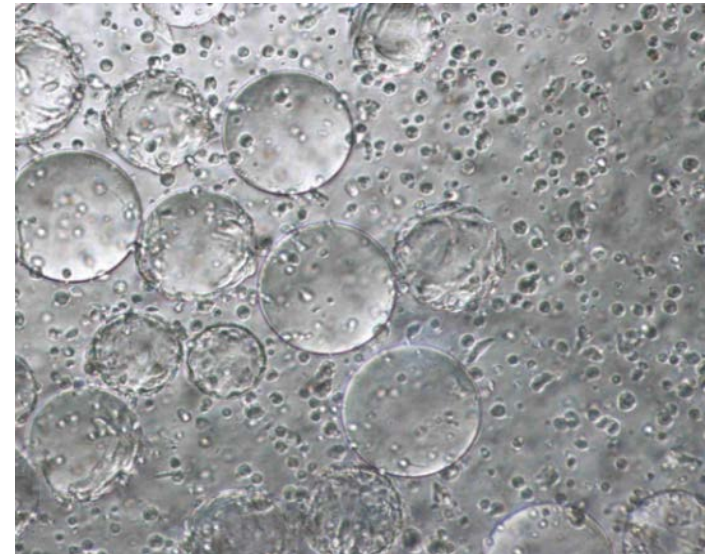


HA concentration at harvest was close to 12 μg/mL
and the virus concentration was $> 10^9$ infective units/mL




Purification Workflow

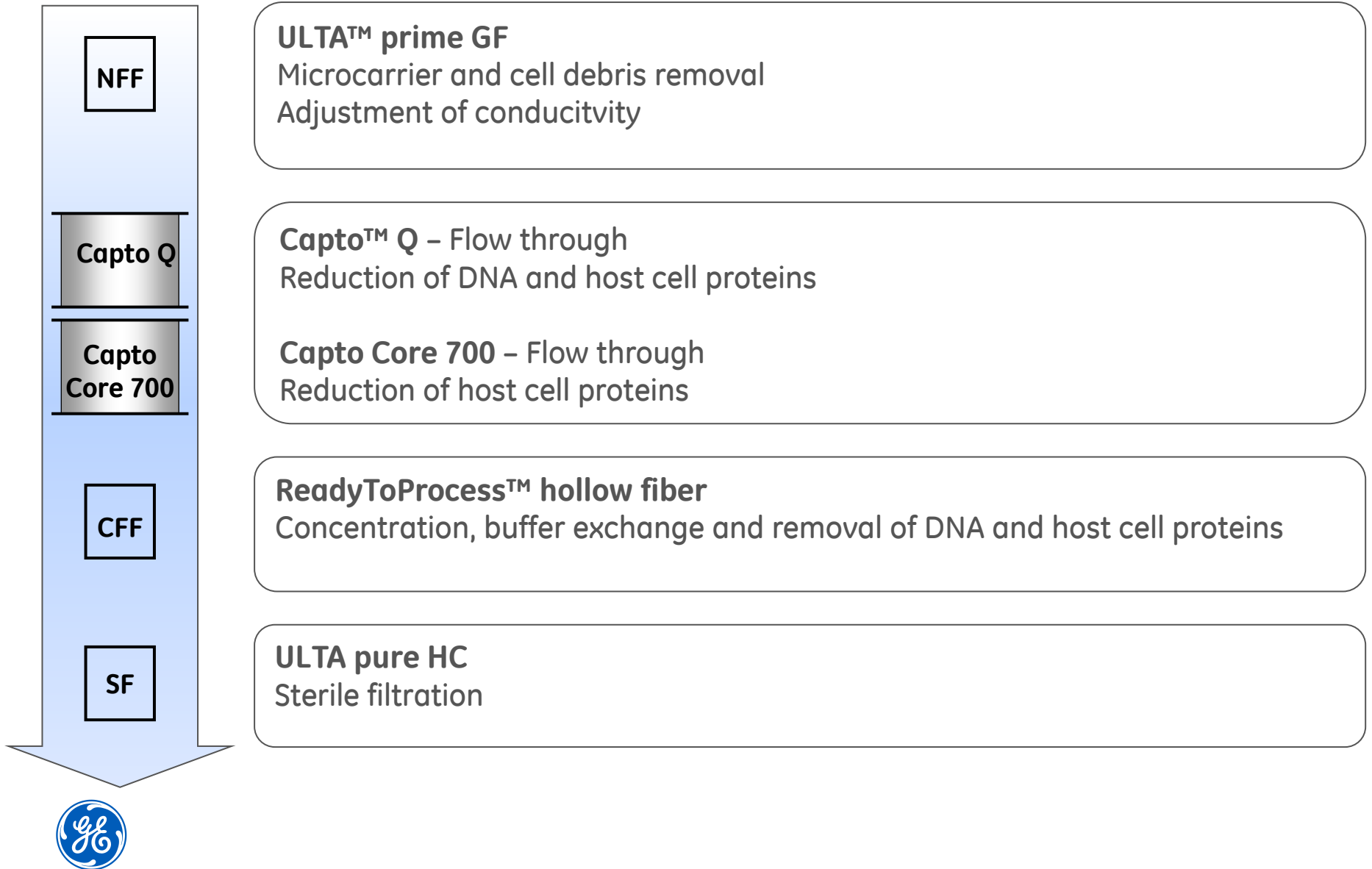


Capto™ Core 700

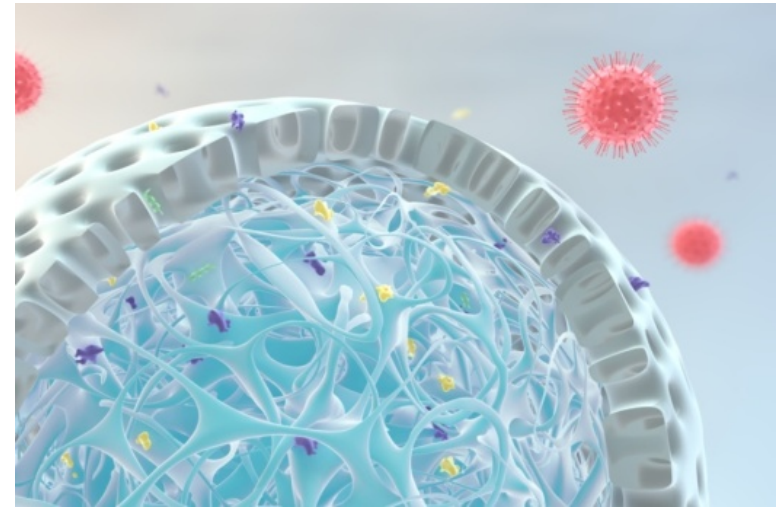

- Shell: Cross-linked native agarose
- Inner core ligand: Octylamine
- Large entites pass outside of the bead while small entities (< 700 kDa) are captured in the core
- Binding not affected by salt , pH or common buffer systems
- Enables flexible unit or serial operation

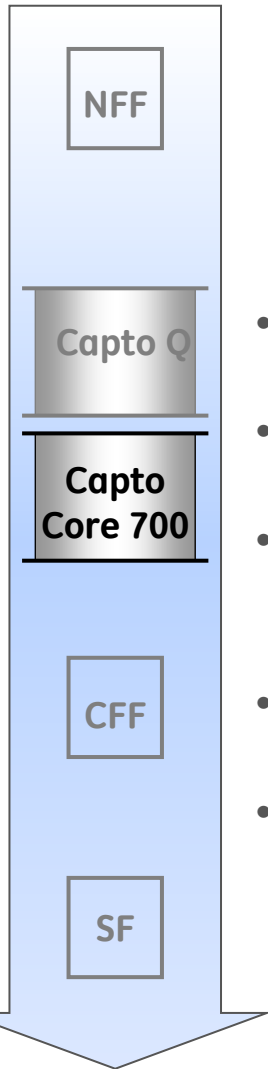


Purification results

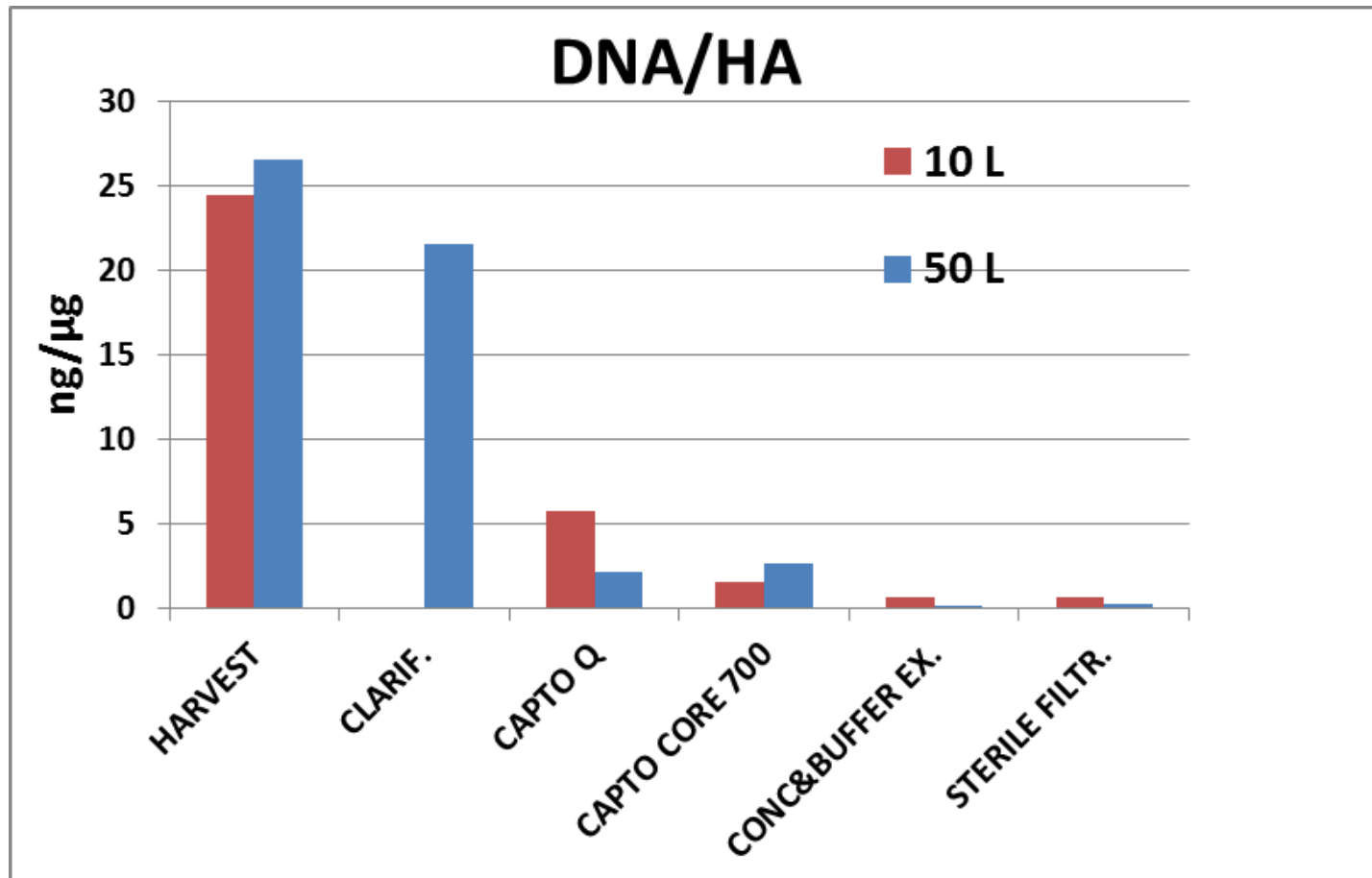

Capto™ Q: Reduces host cell DNA




Purification results

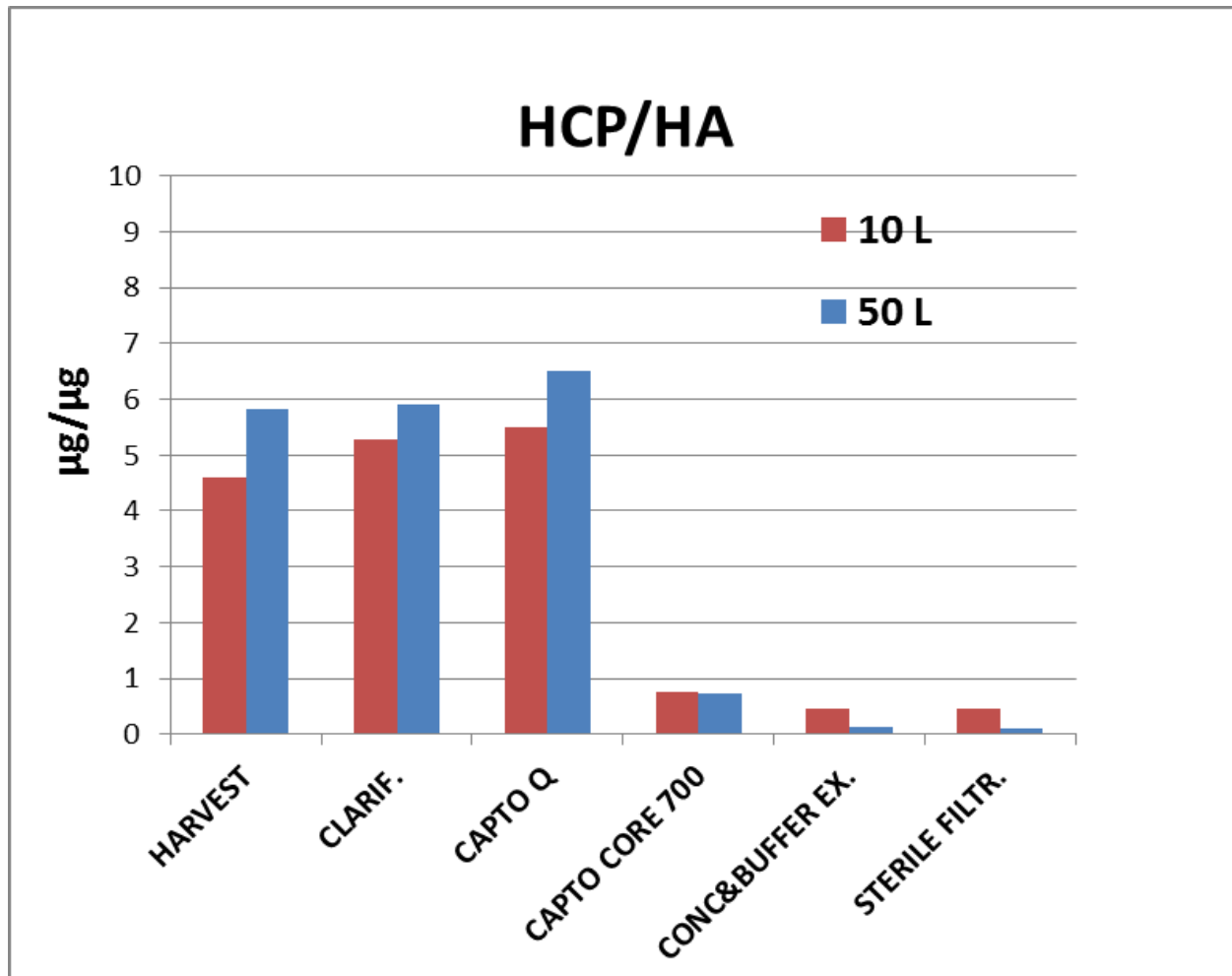

Capto™ Core 700: Reduces host cell protein



Virus infectivity

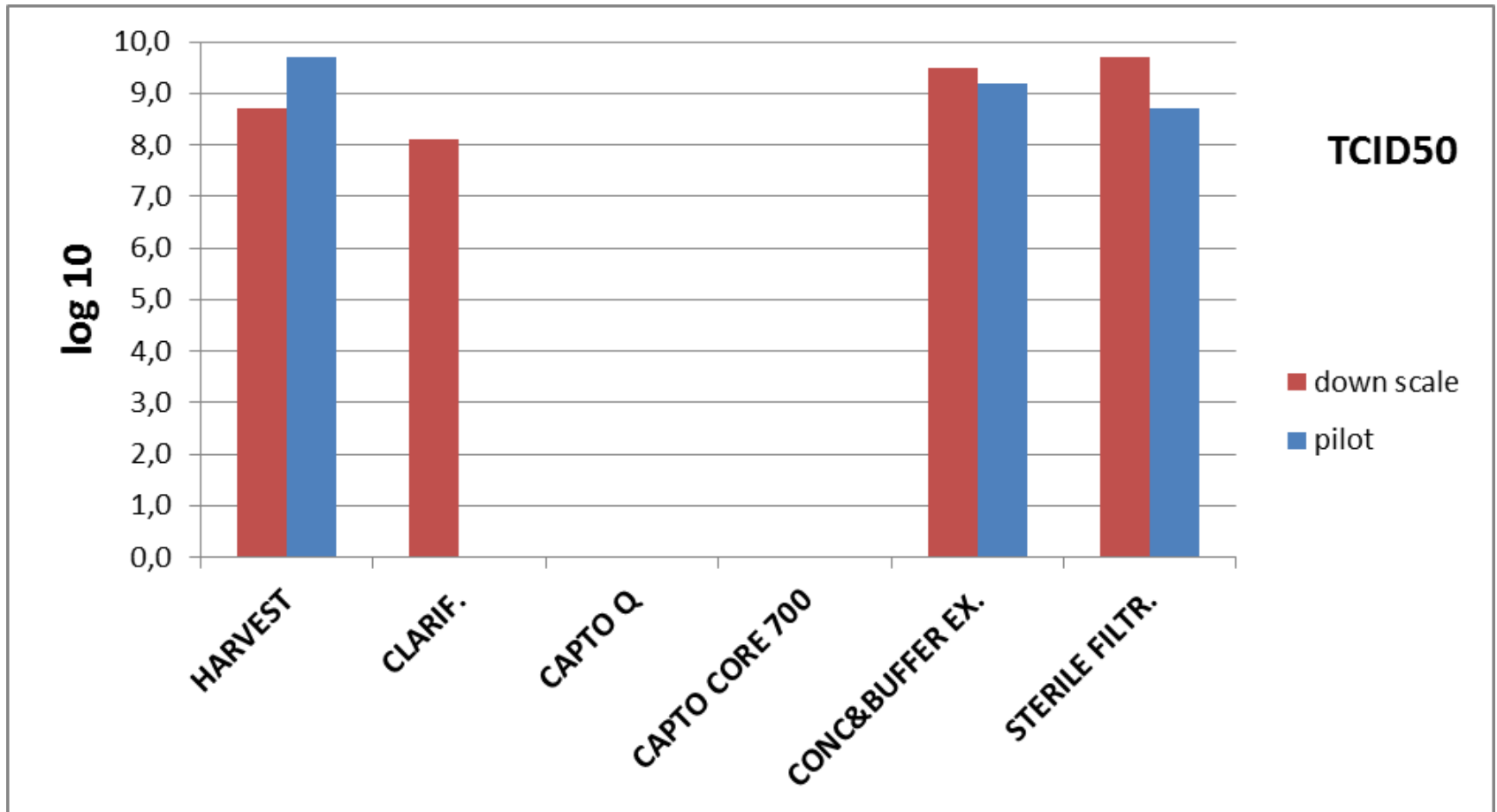

Process does not impair virus infectivity



Process summary

Estimation of doses per liter harvest, compared with WHO guidelines for protein and DNA impurities in influenza vaccine

	Split-inactivated vaccine ¹	Nasal LAIV ²
Scale-up output/L harvest	175 doses á 15 µg HA	3075 doses á 10 ⁷ TCID ₅₀ units
Harvest volume to produce 10 ⁶ doses	5760 L	325 L
Protein impurity ³	30 µg protein/15 µg HA	1.5 µg protein/10 ⁷ TCID ₅₀ units
DNA impurity ⁴	3.0 ng/15 µg HA	0.15 ng/10 ⁷ TCID ₅₀ units

¹ Split-inactivated vaccine contains 3 strains á 15 µg/HA (e.g., 3 × 15 = > 45 µg HA/dose á 0.5 mL).

² Comparison is based on a commercially available specification for a nasal LAIV. A dose of 0.2 mL contains 10⁷ fluorescent focus units, which is assumed to be equal to TCID₅₀ titer.

³ WHO guideline for protein impurity: max. 100 µg protein/strain

⁴ WHO guideline for DNA impurity: < 10 ng DNA/dose = 3.3 ng DNA/15 µg HA.

Assuming a recovery of 25% for the overall process and a dose requirement of 10⁷ TCID₅₀, more than 1.5 million doses of monovalent live attenuated influenza vaccine could be produced from a 50 L cell culture



Yellow fever virus propagation – from eggs to cells



GMP manufacturing of viral vaccine



Xcellerex™ XDR-50 bioreactor

Vero cells (WHO-10-87)


- Cytodex™ 1 microcarrier
- Serum free, animal component-free medium

Yellow fever virus 17D



Virus production drain down refeed

PFU equivalents from Eng and GMP bioreactor runs

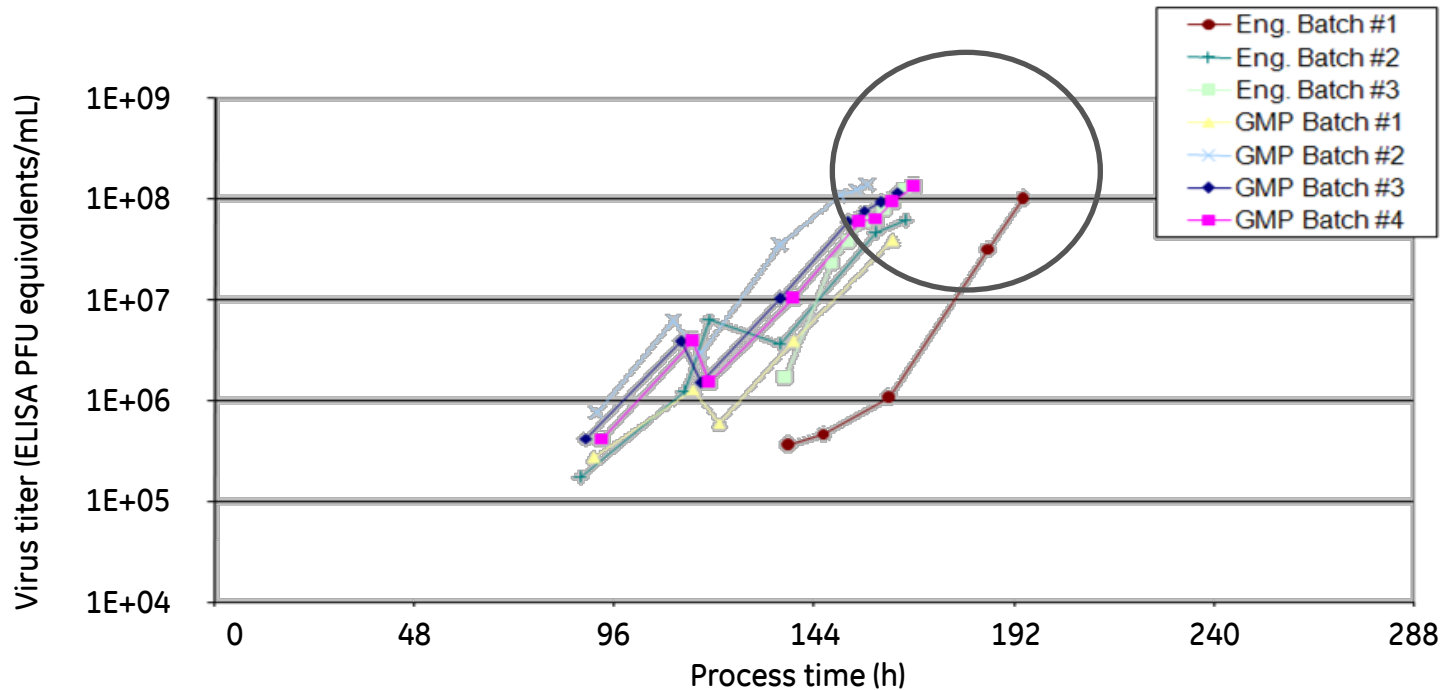


The process consistency was high and virus titers were similar between runs



Virus propagation and release of HCP

HCP content after ELISA analysis

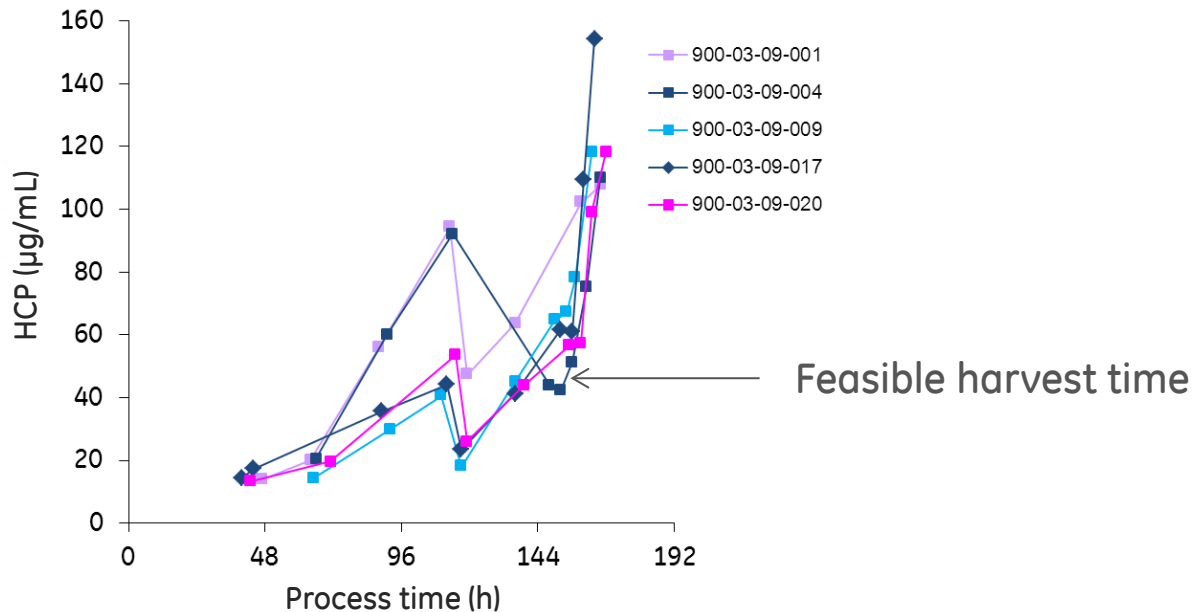


A feasible time for harvest is before the HCP peaks,
to facilitate downstream processing



Conclusions



Conclusions

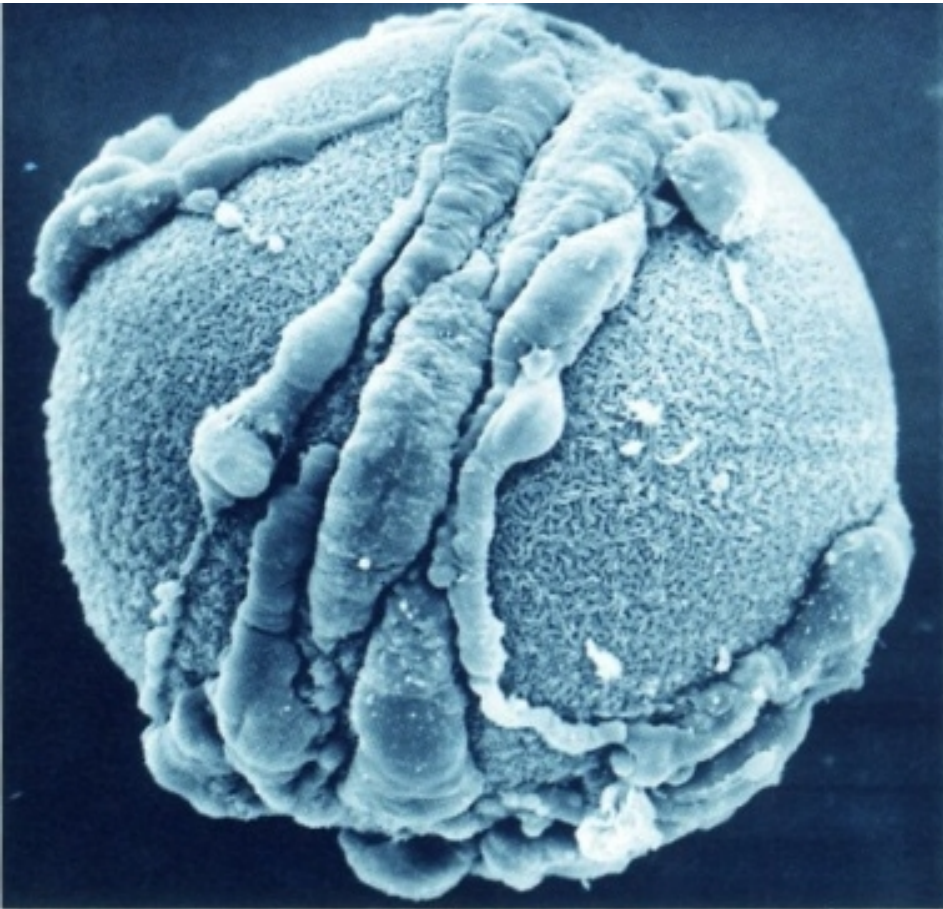

By modernizing legacy vaccine processes there can be improvements in:

Yield

Quality

Scale-up

Cost efficiency

Less hazards



Acknowledgement

Dr. Tian, Yang Michelle (Fast Trak, GE Healthcare Life Science)

Dr. Zheng, Yongjian (Fast Trak, GE Healthcare Life Science)

Dr. Lu, Lifang Lillian (Fast Trak, GE Healthcare Life Science)

Tushar, Joglekar (Fast Trak India, GE Healthcare Life Science)

Chen, Yongshuo (GE Healthcare Life Science)

Liao, Yi Jazmine (GE Healthcare Life Science)

Li Xiaoling (Fast Trak, GE Healthcare Life Science)

Dr. Sui, Lili (Fast Trak, GE Healthcare Life Science)

Luo, Liang (GE Healthcare Life Science)

Xie, Hongyan (GE Healthcare Life Science)

Dr. Jia, Guodong (Fast Trak, GE Healthcare Life Science)

Yang, Hongyan (Fast Trak, GE Healthcare Life Science)

Dr. Lu, Lifang (Fast Trak, GE Healthcare Life Science)

Wang, Jiyuan (Fast Trak, GE Healthcare Life Science)

Yang, Jianjun (Fast Trak, GE Healthcare Life Science)

Du, Qizhen (Fast Trak, GE Healthcare Life Science)

Dr. Lundgren, Mats (GE Healthcare Life Science)

Dr. Berg, Mikael (GE Healthcare Life Science)

Larsson, Anders (GE Healthcare Life Science)



Acknowledgement



Uppsala team:

Lena Sandberg
Peder Bergvall
Camilla Estmer- Nilsson
Christian Kaisermayer
Ann-Christin Magnusson
Christine Sund-Lundström
Therese Lundberg
Mats Lundgren
Jakob Liderfelt
Annika Morrison
Elisabet Linde
Åsa Eriksson
Johanna Tschöp

Marlborough team:

Geoff Hodge
Patrick Guertin
Ed Hayman
Steve Turbayne
Casey Cunningham
Greg Mantenuto
A. J. Theriault
Chris Mach



Thank you!

www.gelifesciences.com

GE, imagination at work, GE monogram, and Global Star are trademarks of General Electric Company.

ÄKTA, Biacore, Capto, Cytodex, HyClone, ReadyCircuit, ReadyToProcess, Sephadex, ULTA , WAVE Bioreactor, and Xcellerex are trademarks of GE Healthcare Company or one of its subsidiaries.

HyperFlask is a trademark of Corning. Pluronic is a trademark of BASF.

© 2013-2015 General Electric Company – All rights reserved. First published Apr. 2013.

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare which supplies them. A copy of these terms and conditions is available on request. Contact your local GE Healthcare representative for the most current information.

GE Healthcare Bio-Sciences AB
Björkgatan 30
751 84 Uppsala
Sweden



