



Centre for Biopharmaceutical Excellence



Centre for Biopharmaceutical Excellence

ICH Q10 Pharmaceutical Quality System (PQS)


© CBE Pty Ltd

This training program is copyright to CBE Pty Ltd and may not be modified, reproduced, sold, loaned, hired or traded in any form without its express written permission.

© CBE – DCVMN 010 V4

Introduction

1



Centre for Biopharmaceutical Excellence

PQS - Module Outcomes

On completion of this module the participant should be able to:

- Interpret ICH Q10 expectations
- Develop a Quality Manual (QM) template and Quality Policies (POL)
- Design a pharmaceutical quality system using the QM template
- Strengthen GMP compliance systems

© CBE – DCVMN 010 V4

Introduction

2

CBE
Centre for Pharmaceutical Excellence

PQS Module Topics

CBE

Basics of Quality and Compliance

CBE

How does a Pharmaceutical Quality System fit together ? – ICH Q 10

CBE

Deviations

CBE

CAPA and Continuous Improvement

© CBE – DCVMN 010 V4 Introduction

CBE
Centre for Pharmaceutical Excellence

Increasing Industry Complexity

COMPLIANCE

© CBE – DCVMN 010 V4 Deviations


Major International Codes of GMP

- **PIC/S Guide to Good Manufacturing Practices - PE 009 - 2014**
- **EU Guide to Good Manufacturing Practices (Eudralex Ch 4)**
- **World Health Organisation (WHO) c GMPs**
- **United States - FDA CFRs Part 21**
 - CFR 210/211 for Drugs and Biologics - current GMPs
 - CFR 820 Quality Systems for Medical Devices - current GMPs
- **ICH Q7 GMP for Active Pharmaceutical Ingredients**
- **Canadian cGMP (aligned with PICs)**
- **ISO 13485 : 2003 - Medical Devices**
- **ICH Guidance Documents – Technical Standards**

Some Useful Reference Documents

- **EU/PICs/TGA cGMPs – Chapter 1 – Quality Management**
- **ICH Q10 - Pharmaceutical Quality System**
- **ICH Q8 – Pharmaceutical Product Development**
- **ICH Q9 - Risk Management in Pharmaceuticals**
- **FDA Quality Systems Approach to Pharmaceutical CGMP Regulations (9/2006)**
- **GHTF - GHTF/SG3/N15R8 - Implementation of risk management within a Quality Management System**

Significant Changes 40 Year History of Pharmaceutical Quality Management

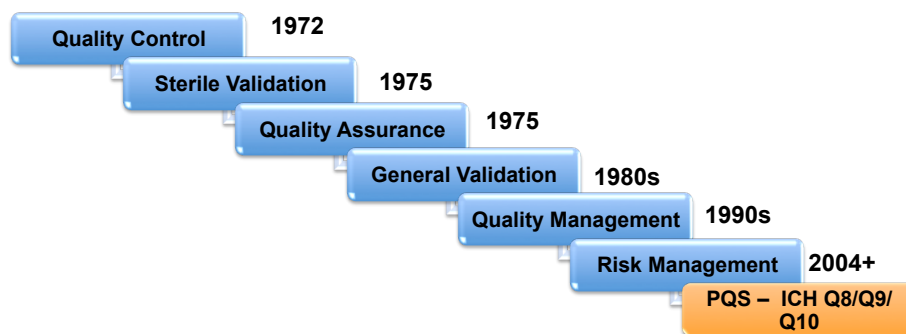


© CBE – DCVMN 010 V4

Introduction

7

Quality Control and Sampling

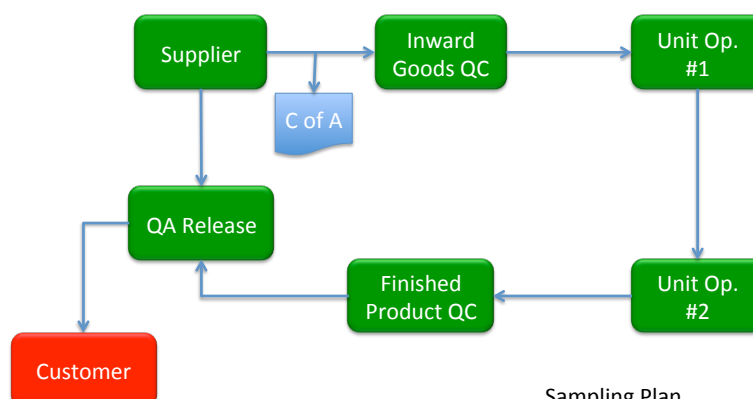

Sampling Plan
Lot Size = 1000
Sample size = 20 (2%)

© CBE – DCVMN 010 V4

Introduction

8

Quality and Compliance

- **Quality refers to:**
 - **Product Quality** – meeting agreed specifications
 - **Quality Systems** – planned and deployed processes (systems) used to monitor, report and take corrective actions
- **GMP Compliance refers to:**
 - Identification, documentation and deployment of GMP obligations
 - Ongoing verification that GMP obligations are being met, or not.

An Overview of ICH Q10 Pharmaceutical Quality System

ICH Q10 - Pharmaceutical Quality System

- Based on ISO 9000/ISO13485/CFR 820 systems model
- Compliments ICH Q8 and ICH Q9
- Applies across the product life-cycle
- Consistent with GMPs - not intended to add new expectations to regulations and compliance
- Applies to APIs, drug products and biotechnology
- Strengthens the link between product development and manufacturing activities

Pharmaceutical Quality System, Quality Assurance, GMP and Quality Control

Pharmaceutical Quality System (ICH Q10)



ICH Q10 - Pharmaceutical Quality System

CBE
Center for Pharmaceutical Excellence




13



A Hierarchy QMS/GMP Documentation

CBE
Center for Pharmaceutical Excellence



© CBE - D

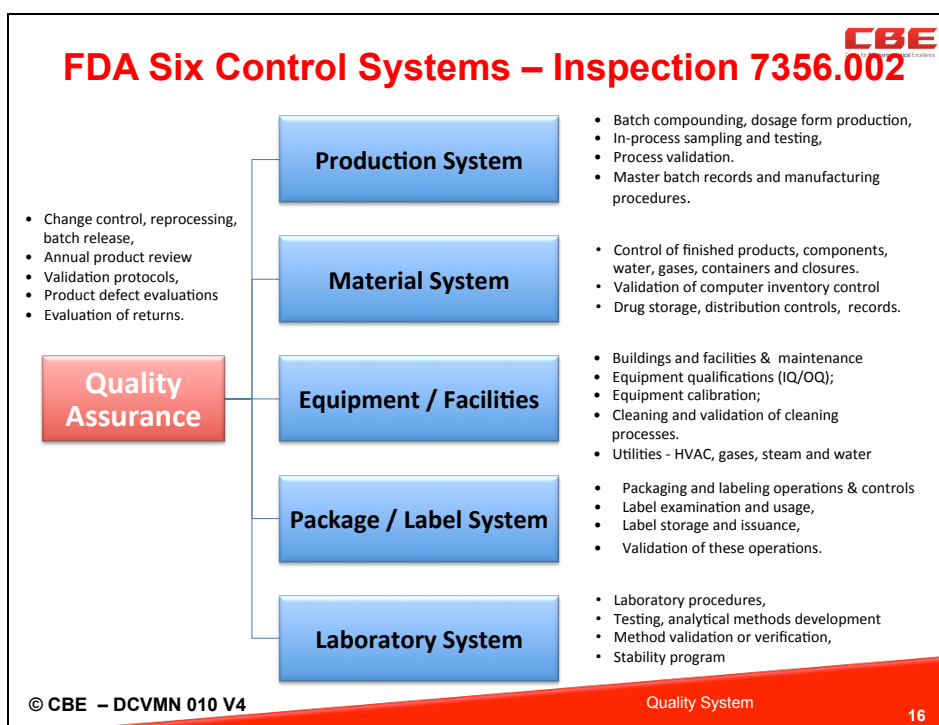

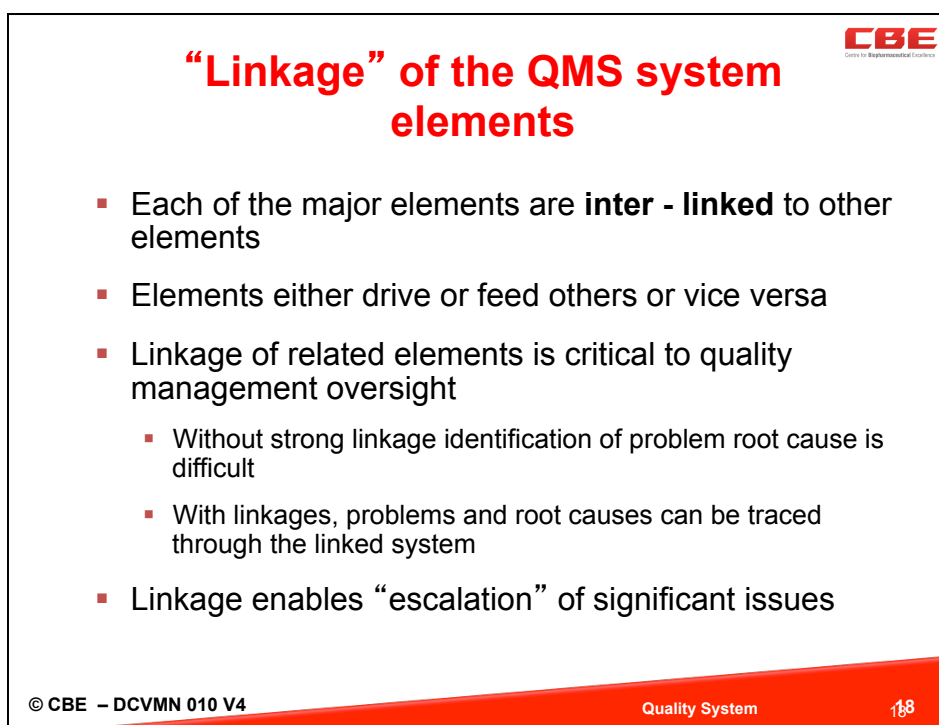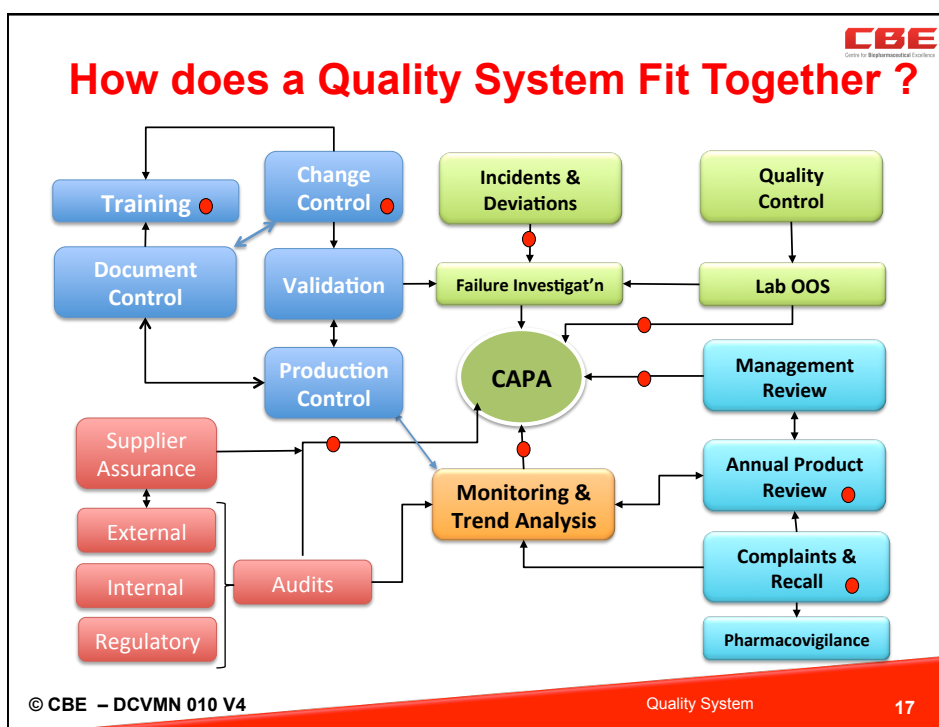


ICH Q10 - Some Important Principles



- The size and complexity of the company's activities should be taken into consideration when developing a new pharmaceutical quality system or modifying an existing one.
- While some aspects of the pharmaceutical quality system can be company-wide and others site-specific, the effectiveness of the implementation of the pharmaceutical quality system is normally demonstrated at the site level.

Integration of PQS and GMP Elements in the Quality System

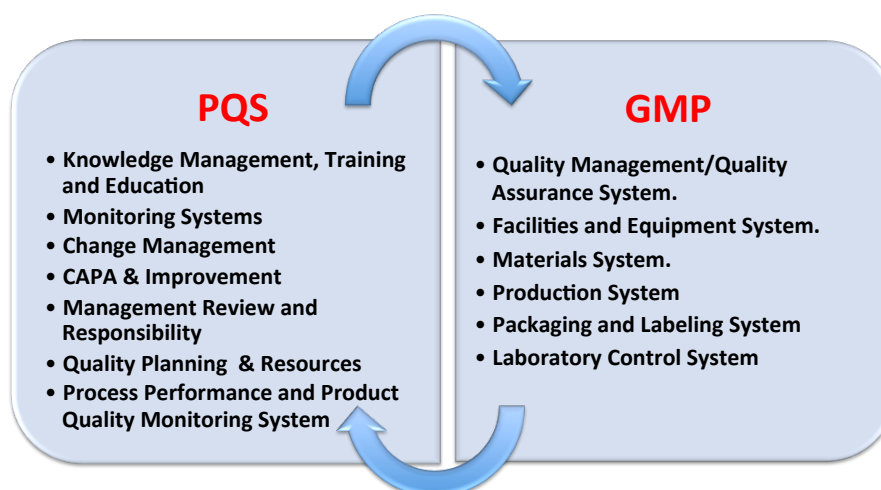

PQS Enablers Quality Risk Management (QRM)



- **Quality risk management**, in line with ICH Q9, provides an essential component of the Quality System.
- QRM enables both effective and efficient practices.
- Application of QRM ensures the quality system is efficient.
- Provides a systematic approach to escalating and prioritising significant events



PQS Enablers - Knowledge Management



- **Knowledge management** means the systematic accumulation of information concerning products so that this knowledge can be leveraged in the future.
- **Knowledge** can be stored in systems such as:
 - Quality Records, including testing, stability studies and reports
 - Registration Dossiers
 - Contracts and Technical Agreements
 - Validation Protocol and Reports
 - Marketplace Events (Complaints, Recalls, Adverse Events etc.)
 - Annual Product Quality Reviews (PQRs)
 - Process control, significant deviations and changes.



ICH Q10 – Pharmaceutical Quality Manual



- A Quality Manual or equivalent documentation approach should be established and should contain the description of the pharmaceutical quality system.
- The description should include:
 - i) The quality policy
 - ii) The scope of the pharmaceutical quality system.
 - iii) Identification of the **processes** within the pharmaceutical quality system, as well as their **sequences, linkages** and **inter-dependencies**.
- Process maps and flow charts can be useful tools to facilitate depicting these in a visual manner.

Example Chapters of a Quality Manual



- MANAGEMENT REVIEW AND RESPONSIBILITY
- DESCRIPTION OF THE QUALITY SYSTEM
- QUALITY PLANNING AND RESOURCE MANAGEMENT
- TRAINING AND EDUCATION
- PRODUCT DEVELOPMENT AND PRODUCT REGISTRATION
- QUALITY ASSURANCE AND COMPLIANCE PROGRAMS
- MONITORING PROGRAMS
 - CAPA / QUALITY AUDITS / PRODUCT QUALITY REVIEWS
 - MARKETPLACE MONITORING: COMPLAINTS AND PHARMACOVIGILANCE PROGRAMS
- QUALITY RISK MANAGEMENT
- KNOWLEDGE MANAGEMENT



Example Chapters of a Quality Manual

- **MATERIAL CONTROL SYSTEM:**
 - SUPPLY CHAIN INTEGRITY AND SUPPLIER ASSURANCE
- **PRODUCTION SYSTEM**
 - PRODUCT DEVELOPMENT, TECHNOLOGY TRANSFER AND MASTER INSTRUCTIONS
 - PROCESS PERFORMANCE AND PRODUCT QUALITY MONITORING SYSTEM
 - PROCESS VALIDATION
 - PACKAGING AND LABELLING
 - DEVIATIONS, INVESTIGATIONS AND NON-CONFORMING PRODUCT
- **LABORATORY CONTROL SYSTEM**
- **CHANGE MANAGEMENT / VALIDATION**
- **FACILITIES AND EQUIPMENT SYSTEM / CLEANING/ CONTAMINATION CONTROL AND COMPUTERISED SYSTEMS**



PQS Objectives

- Set objectives and clear performance metrics;
 - Objectives are reviewed annually
 - Metrics are measured, time-related and indicate the level of performance required;
 - Metrics are reviewed regularly;
- Define how compliance obligations are embedded in operational practices and procedures.
- Address processes for identifying, reporting and responding to compliance failures.

FDA View on Quality Metrics

Indicator	Metric
Lot Acceptance Rate	Number of lots rejected in a year / number of lots produced
Right First Time Rate	Number of deviations / lot
Complaint Rate	Number valid complaints/number of lots released per year
Invalidated (OOS) Rate	Number of OOS test results invalidated /tests performed
Annual Product Review (APR) on Time Rate	Number of APRs generated within 30 days of annual due date
Management Engagement	Most senior manager that signed each annual product review
Process capability or performance index	Whether performed for each critical quality attribute as part of that product's APR.
Corrective and Preventative Action (CAPA) Rate	Number of CAPAs that were initiated due to an APR, divided by the total number of APRs generated.

ICH Q10 – Management Review

- Senior management should be responsible for pharmaceutical quality system governance through management review to ensure its continuing suitability and effectiveness.
- Management should assess the conclusions of periodic reviews of process performance and product quality and of the pharmaceutical quality system.

Commitment to Quality by Management

- Governing Body and CEO are engaged
- The Quality Policy is aligned with business objectives
- Compliance and Quality obligations are embedded in position responsibilities
- Resources are allocated to Quality / Compliance
- Top level engagement in compliance/quality metrics and reviews


PQS, Compliance and Quality Organisation

- Assign compliance and quality responsibilities to individual managers – set out in position descriptions
- Ensure all management “walk the talk”
- Appoint a senior Compliance/ Quality executive:
 - direct access to the Board/CEO
 - Access to expert advice (internal and external)
 - Establish compliance/quality objectives and KPIs
- Ensure compliance/quality function has the authority to act

ICH Q10 - Management Reviews – should include

- A timely and effective **escalation** process to senior management;
- Measures of customer satisfaction - complaints and recalls;
- Conclusions of process performance and product quality monitoring;
- The effectiveness of process and product changes including those arising from CAPA;
- Any follow-up actions from previous reviews;

Management Reviews, Trend Analysis and Feedback

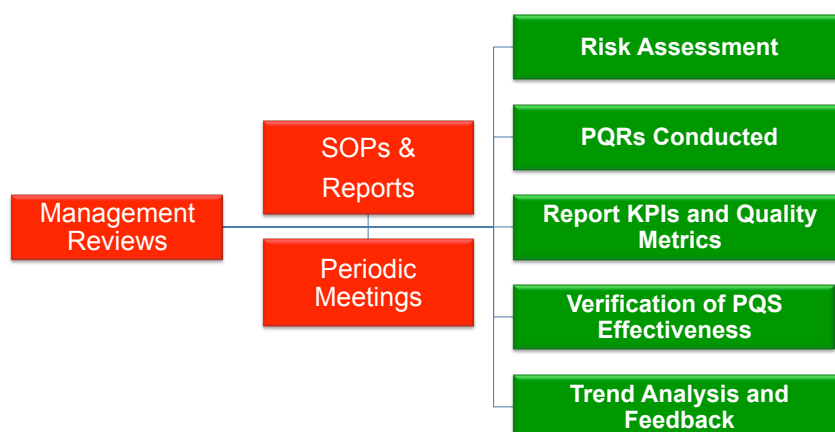

ICH Q10 Quality System

Continual Improvement of Process Performance & Product Quality



- **Process performance and product quality monitoring system:**
 - Well defined systems
 - Process control
 - Identification of improvement areas
- **Corrective action and preventive action (CAPA) system**
 - In place and effectiveness evaluated
 - Focus on Continuous Improvement
- **Change management system:**
 - QA oversight
 - Utilizes science and risk-based assessment
- **Management review of process performance and product quality**
 - Periodic reviews of performance against metrics
 - Supports continual improvement

© CBE – DCVMN 010 V4

Quality System

ICH Q10 Quality System:

Continual Improvement



- **Management Review of the Pharmaceutical Quality System**
 - Measurement of achievement of QS objectives
 - Assessment of Metrics
- **Monitoring of Internal and External Factors impacting the QS**
 - Emerging regulations, guidance and quality issues
 - Innovations
 - Changes in business strategies and objectives.
- **Outcomes of Management Review and Monitoring**
 - (Re)allocation of resources and/or personnel training
 - Timely and effective communication of the results

© CBE – DCVMN 010 V4

Quality System

34

Compliance and Improvement




Element	Compliance is a cost \$	Regulation Driven	Improved Compliance	Integrated Compliance	Competitive Advantage
Quality Focus	Testing (QC)	GMPs	Processes and Systems	ICH Q8, 9, 10 Started	ICH Q8,9,10 Embedded
CAPA	Correction (Reactive)	Corrective Action	Prevent. Action (RCA)	Management Reviews	Drive down COQ
Continuous Improvement	Absent	Event Driven (Reactive)	QA Focus (Predictive)	Operations Focus	Company Wide - Part of Culture
Compliance	QA/QCs role - minimal Audits	Compliance / GMP Audits	PQS Systems driven audits	Prepare for Regulatory Audits	Welcome External Feedback
Knowledge & Training	Basic GMP Training	Knowledge is anecdotal	Systematic Training Evaluated	Knowledge is Documented & Organised	Knowledge is Leveraged

© CBE – DCVMN 010 V4

35

Deviations, Investigations and CAPA

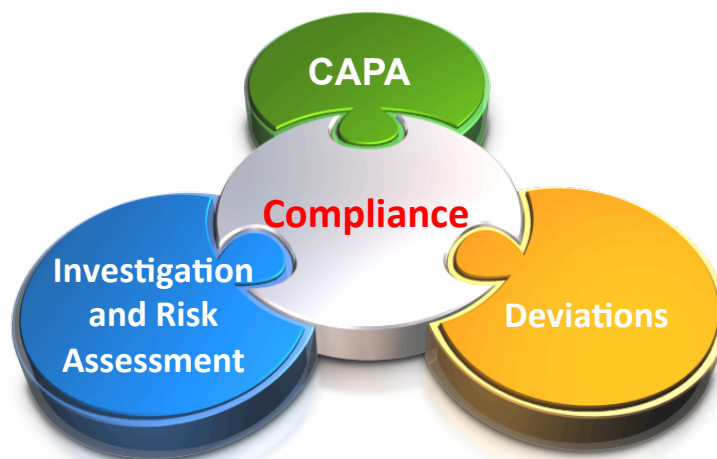

© CBE – DCVMN 010 V4

Deviations

36

Top Critical/Major Defect Areas

Most Frequent Defect Categories Observed		
Rank	Defect Category	Percentage of Critical / Major Deficiencies with this Defect Category
1	Investigation of anomalies	6.5%
2	Quality management	5.5%
3	Investigation of anomalies – CAPA	4.7%
=4	Contamination, chemical/physical (or potential for)	3.7%
=4	Supplier and contractor audit	3.7%
6	Quality management – change control	3.6%
7	Documentation - procedures/PSF/TAs	2.7%
7	Personnel issues – training	2.7%
=9	Design and maintenance of equipment	2.6%
=9	Documentation – manufacturing	2.6%
=9	Finished product testing - chemical	2.6%

The Defect Categories presented above account for 40.9% of all Critical / Major Deficiencies raised in 2013.

© CBE – DCVMN 010 V4

37

Manufacturers' Obligations

- **Oversight:** If a deviation occur, it should be approved in writing by a competent person, with the involvement of the Quality Control Department when appropriate;
- **Release:** Deviations should be resolved before release
- **Stability:** Significant batch deviations may invoke a stability trial
- **PQR:** A review of all significant deviations, their related investigations, and the effectiveness of resultant CAPA taken.

© CBE – DCVMN 010 V4



Deviations

PICs cGMP Expectations

- Any significant deviations are fully recorded and investigated; (GMP 1.3 iv.)
- Product assessment an assessment of deviations from specified procedures; (GMP 1.3 vi.)
- Notes on special problems including details, with signed authorisation for any deviation from the Manufacturing Formula, Processing Instructions and Packaging Instructions (GMP 4.17 (i.)/ GMP 4.18 (h))
- If a deviation occur, it should be approved in writing by a competent person, with the involvement of the Quality Control Department when appropriate. (GMP 5.15)

PICs cGMP Expectations

- Any significant deviation from the expected yield should be recorded and investigated. (GMP 5.39)
- an on-going stability study should be conducted after any significant change or significant deviation to the process or package. (GMP 6.30)
- The Competent Authorities should be informed if a manufacturer is considering action following possibly faulty manufacture, product deterioration, detection of counterfeiting or any other serious quality problems with a product. (GMP 8.8)

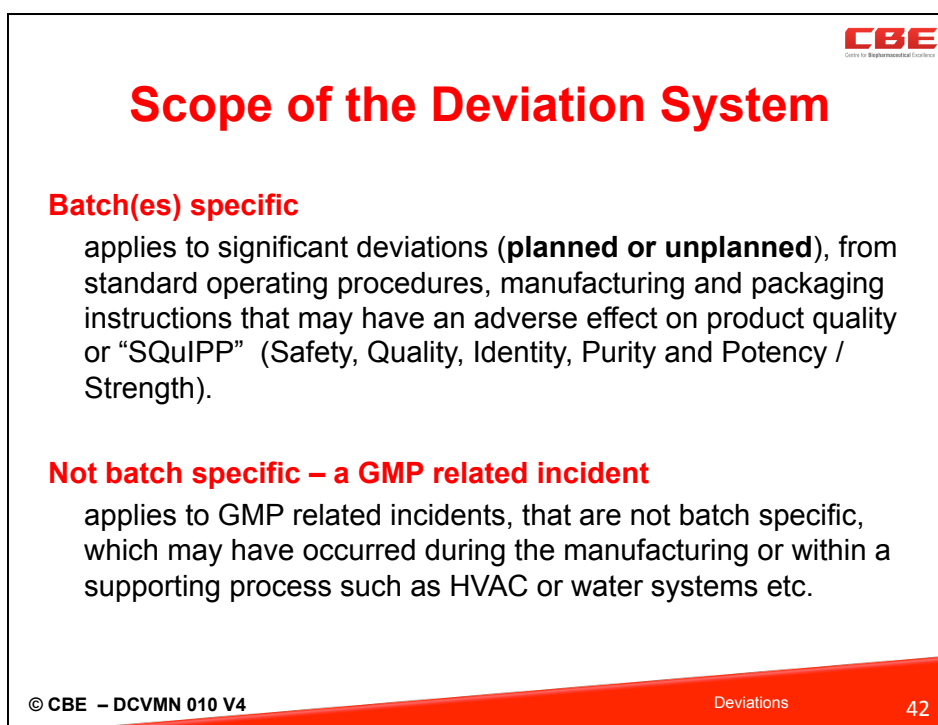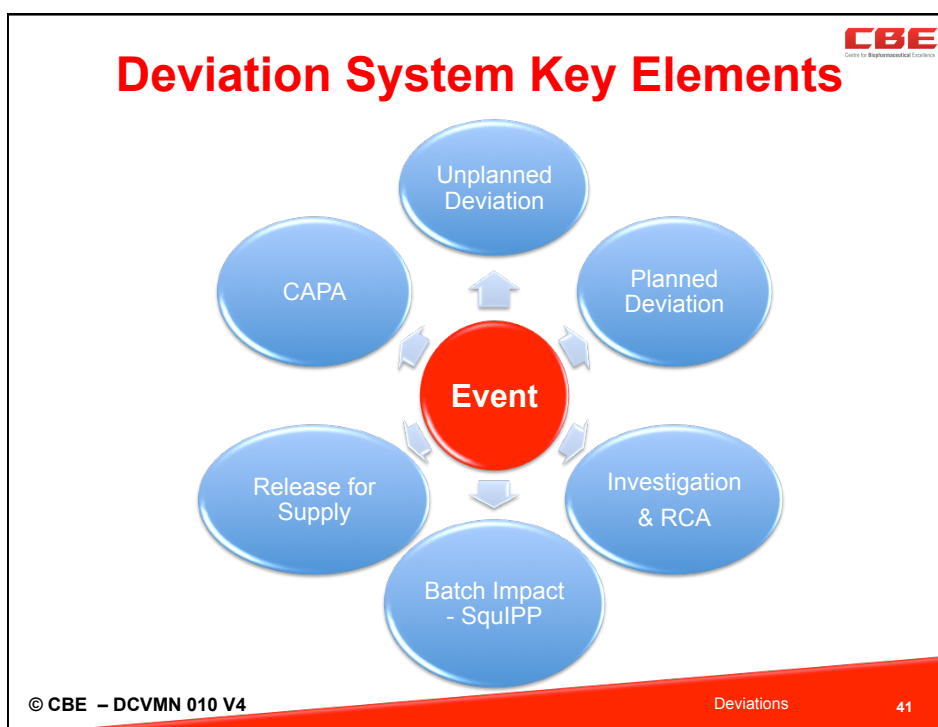

Deviation System also does apply to



- Maintenance and calibration – relating to GMP equipment and services
- Confirmed Out of Specification (OOS) events
- Laboratory procedures and test methods
- Stability failures
- Environmental monitoring and other GMP service excursions from action limits
- Supply chain / raw materials integrity
- Concurrent process validations and cleaning validations
- Phase III clinical trials material manufacture

Deviation System does not apply to



- Audit observations
- Product complaints and adverse events
- Returns and recalls
- Prospective qualifications and validations (these are handled within the Validation Master Plan procedures)
- Clinical Trial materials - Phase I / II manufacture

Healthy Deviation Management Environment



- Staff feel able to raise a deviation without blame
- Deviations are expected – its how we manage them that counts
- Good communication and judgment around when to report, or not – seek advice
- Constructive use of investigation and risk assessment tools



Alert and Responsive to Events

Responsibility of QA




- Approval of planned deviations **before their implementation**;
- Classification of the deviation on the basis of Risk
- Overseeing a deviation investigation and review of any investigation / impact assessment report
- Filing completed deviation and incident reports.
- Deciding if a CAPA is required, or not
- Assessing subsequent corrective actions and investigation details
- Reviewing a deviation or incident report at point of release for use or for supply;
- Disposition of the product or material
- Updating and maintaining the Deviation/Event register

Deviation Decisions

- Should all Quality related “Events” be recorded ?
- Should all Events be referred to QA ?
- When does an event become a GMP deviation ?
- How is a “Significant” deviation defined ?
- Should all deviations be investigated ?
- How do we know its significant if its not investigated ?
- Should all investigations be documented / risk assessed ?
- Should CAPA be applied to all investigation outcomes

Filtration and Escalation Approach

- Start with a GMP related “Quality Event”
- Record on an Event Log
- Assess its potential significance
- Escalate to Deviation, or not
- Commence investigation

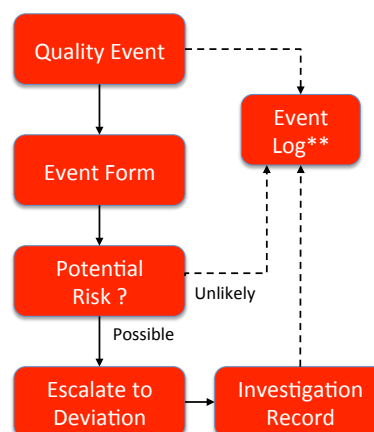

** Periodic and Annual Review

Categorise Event Categories for Trending



- ☐ Excursion from MBR
- ☐ Excursion from SOP
- ☐ Excursion from Test Method
- ☐ EM Excursion
- ☐ Equipment Breakdown
- ☐ Facility Breakdown
- ☐ Materials / Components
- ☐ Other

Example Only

© CBE – DCVMN 010 V4

Deviations

Checklist to Prompt Preliminary Risk Decision (Potential Risk ?)



1	Likely the event could impact SQIIPP ? (Safety, Quality, Identity, Purity, Potency)	Yes No Unsure ?
2	Does the event result in an excursion from registered details for this product ?	Yes No Unsure ?
3	Likely the event could cause physical contamination or cross contamination ?	Yes No Unsure ?
4	Likely the event could cause loss of identity or traceability ?	Yes No Unsure ?
5	Likely the event could result in an out of specification result, if tested ?	Yes No Unsure ?
6	Likely the event could affect product quality or stability in the marketplace ?	Yes No Unsure ?
7	Is the event related to a GMP non-conformance or outside the "validated state" ?	Yes No Unsure ?
8	Likely the event has compromised a CPP or a CQA ?	Yes No Unsure ?

© CBE – DCVMN 010 V4

Recording and Evaluating Deviations

Most important to record the deviation quickly and accurately.

Record

- Date / time/ process step and stage of processing (pallet #)
- Batch #(s) and Item #(s)
- Equipment, process line and operator(s)
- Sequence of events causing the deviation
- How the deviation was identified
- What immediate action was taken (Containment)

Evaluation is very dependent on good records

- Line and product trend history
- Manufacturing batch records and line logs
- Level of in-process controls

Report Number:	QE/DR -	Title
Refer to SOP 018030 when completing this form.		If this form is typed initial here to confirm version control has been checked:
1. Quality Event Notice and Assessment <input type="checkbox"/> Planned Deviation** <input type="checkbox"/> Unplanned Deviation		
Origin (✓)	Description of the Event	
<input type="checkbox"/> Excursion from MBR <input type="checkbox"/> Excursion from SOP <input type="checkbox"/> Excursion from Test Method <input type="checkbox"/> EM Excursion <input type="checkbox"/> Equipment Breakdown <input type="checkbox"/> Facility Breakdown <input type="checkbox"/> Materials / Components <input type="checkbox"/> Other	(Batch # / Date / Process Line, Personnel etc.) or description of planned deviation and associated control actions.	
<input type="checkbox"/> Not Significant <input type="checkbox"/> Potentially Significant (Progress to deviation)	Reasons for Conclusions:	
** Planned deviations require the prior approval of the QA Manager		
Signed / Dated by QA Manager		

Record the event in the Quality Event Trend Register – If the event has been assessed as Not Significant record as QE or, if assessed as Potentially Significant then record as a Deviation (DEV) and continue with the completion of this form.

Report Number:	CAPA -	Title					
Refer to SOP 021 when completing this form.		If this form is typed initial here to confirm version control has been checked:					
This record combines (a) failure investigation (b) final risk assessment and (c) CAPA for any significant quality related event. It may be used stand alone or as a result of preliminary review of deviations, audits, complaints, vendor assessments or laboratory OOS.							
Short Description of the Issue		Reference #: **					
** Deviation #, Internal audit #, Complaint # or external audit # etc.							
<table border="1"> <tr> <th colspan="2">Origin (✓)</th> </tr> <tr> <td> <input type="checkbox"/> Unplanned Batch Deviation <input type="checkbox"/> Internal Audit Significant Observation <input type="checkbox"/> External Audit Significant Observation <input type="checkbox"/> Vendor Non-Conformance </td> <td> <input type="checkbox"/> Planned Deviation <input type="checkbox"/> Customer Complaint or Adverse Event <input type="checkbox"/> Laboratory OOS <input type="checkbox"/> Other: </td> </tr> </table>				Origin (✓)		<input type="checkbox"/> Unplanned Batch Deviation <input type="checkbox"/> Internal Audit Significant Observation <input type="checkbox"/> External Audit Significant Observation <input type="checkbox"/> Vendor Non-Conformance	<input type="checkbox"/> Planned Deviation <input type="checkbox"/> Customer Complaint or Adverse Event <input type="checkbox"/> Laboratory OOS <input type="checkbox"/> Other:
Origin (✓)							
<input type="checkbox"/> Unplanned Batch Deviation <input type="checkbox"/> Internal Audit Significant Observation <input type="checkbox"/> External Audit Significant Observation <input type="checkbox"/> Vendor Non-Conformance	<input type="checkbox"/> Planned Deviation <input type="checkbox"/> Customer Complaint or Adverse Event <input type="checkbox"/> Laboratory OOS <input type="checkbox"/> Other:						
If the event impacts product, write the product name, code and batch number here (or attach list):							
Product Name:		Product Code:	Batch No:				
QA Manager (sign)		Date:					
1. Investigation and Test Plans Conduct a formal investigation of the issue. This investigation must be overseen by the QA Manager. Describe the overall investigation plan and any sampling and test/inspection plan.							
Investigation Plan (describe the investigation actions required in things that must be reviewed including any trends and records)		Testing / Inspection Plan (describe any sampling and test plan)					

© CBE - 53


Investigations

- SOP on “Investigations, Root Cause Analysis and CAPA”
- SOP serves all investigations not just Deviations/Events
- Investigation plans should be documented on proformas
 - Manufacturing plan
 - Quality Control / Stability plan
- Plans should examine beyond the immediate event – refer to Event logs and other batches/products
- Assign investigation leads – does not have to be QA

Investigation Tips and Tools

- Not all problems need RCA, or they can be solved simply
- Examine the “scene of the crime”
- Involve an SME
- 7 Management Tools, then 7 Statistical Tools
- Tools should be quickly accessible to users
 - 5 whys / brainstorming
 - Root cause mapping / C&E Diagrams
 - Pareto, Kepner Tregoe, DMAIC
- Last resort – FMEA level approach

Investigation Flowchart

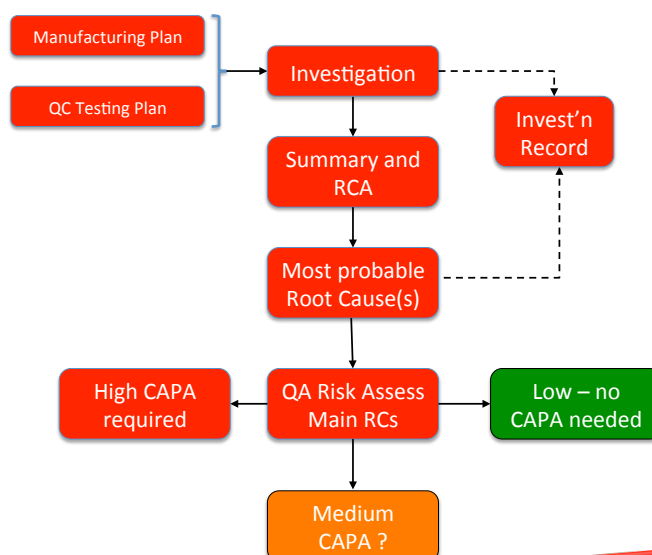

Investigation Tips

- Do it quickly
- Interview Operators
- Root Cause(s) – silver bullet ?
 - Ineffective training
 - Human Error
 - Re-write the SOP
- Trigger Event (generally obvious)
- Underlying Condition (often obscure)



Root Cause Analysis /Investigations – some tips

- Investigate “in the moment”, not with hindsight.
- Be systematic and objective – don’t focus on silver bullet
- Consider “**Look-back**” and “**Look-forward**”
- “Operator Error - Re-train the operator.”
 - Operator error has at least 7 different causes.
 - In a training system that was possibly flawed, to an SOP that may have generated the error ?

Hard Questions in Investigations



- Natural tendency to limit investigations to the batch in question.
- **“Look – back” or “Look – forward”**.
 - Look back – previous batches / products affected
 - Look Forward – likely to repeat the problem in the future – what’s changed ?
- Regulators rightly expect that these potential consequential issues are assessed and documented.
- Not addressing consequential issues is a surefire way to generate a Warning Letter by FDA and criticism from TGA / PICs Inspectors.

Examples



- 1,500 Litres of Vaccine down the drain
- “This batch has glass in it – it shouldn’t be released”
- OOS low potency for biological – repeat the test
- Blow moulded bottle – base uneven

Expanded Investigation (Look back and Look forward)

- When investigating a deviations it is not enough to simply review the event is isolation. The investigations should:
- **Look- back** on past batches that may have been compromised by the deviation under review. Examine batch records, test records other deviation records and complaint records. Look back should determine whether any quarantine or hold on related batches should occur or whether batches released to market should be recalled.
- **Look forward** to try and identify whether future batches may be compromised if no CAPA action is taken. This will determine when processing may recommence and what additional controls may be needed.

Deviation Resolution and Release

- Release: Deviations should be resolved before release of materials or product.
- Does this also require implementation of CAPA ?
 - Correction – Yes
 - CAPA – if possible but not always feasible
- Two point close out for Deviation / CAPA
 - Deviation Closed
 - CAPA Completed

Timeframes for Processing Deviations and Investigations



- Deviations reports should be raised within **2 - 3 working days** of the event occurring and submitted to Quality Assurance.
- Batch/ SQulPP related deviations/incidents must be closed out before any implicated batch is released.
- Close out means that a **batch correction** must be implemented, where warranted.
- All other (non-SqulPP) deviations/incidents should be closed out **within 30 calendar days**.

Assessing Deviation Significance is related to CPPs, CQAs and CSMs



Critical Process Parameter (CPP)

A process parameter whose variability has an impact on a critical quality attribute and therefore should be monitored or controlled to ensure the process produces the desired quality

Critical Quality Attribute (CQA)

A physical, chemical, biological or microbiological property or characteristic that should be within an appropriate limit, range, or distribution to ensure the desired product quality.

Critical Starting Material (CSM)

Critical Quality Attribute(s) of a Starting Material

Outcomes of Investigations

- **Clear SQulPP Impact** – an event or deviation that is likely to have an actual adverse effect on product quality, safety, purity, identity or potency. The deviation is most likely to have an impact on a CPP and/or a CQA.
- **Possible/Probable SQulPP Impact** – an isolated event or deviation from an approved procedure that may have an unknown effect on a product. The deviations may or may not have an impact on a CPP, but is unlikely to have any impact on a CQA.
- **Clear no SQulPP Impact** – an event or deviation that has no actual or a potential adverse effect on product quality, safety or efficacy. The deviation is likely to have no impact on a CPP and/or a CQA.
- **Other** – a deviation from GMP or from a procedure that has very low to no potential impact on product quality or a product CQA / CPP).

Two Things to Keep in Perspective

- The holder of a manufacturing authorisation must manufacture medicinal products so as to ensure that they are fit for their intended use, comply with the requirements of the Marketing Authorisation and do not place patients at risk due to inadequate safety, quality or efficacy.
- The level of effort, formality and documentation of the quality risk management process is commensurate with the level of risk.


CAPA and Regulatory Guidance

- Drug GMPs have been backward looking. Pharma Industry would do well to study the Medical Device requirements.
- PIC/s cGMPs are pretty light on in terms of CAPA expectations – inspectors are not.
- ICH Q10 provides a significant step up in expectation but not mandated yet.
- FDA regularly reference lack of effective CAPA in warning letters.



ICH Q10 - Corrective and Preventive Action

- Should have a system for implementing CAPAs resulting from investigations of:
 - Complaints and Recalls
 - Product Rejections and Non-conformances
 - Deviations
 - Audits & Regulatory inspection findings
 - Trends from process performance and product quality monitoring
- **“The level of effort and formality of investigation depends on the level of risk”**

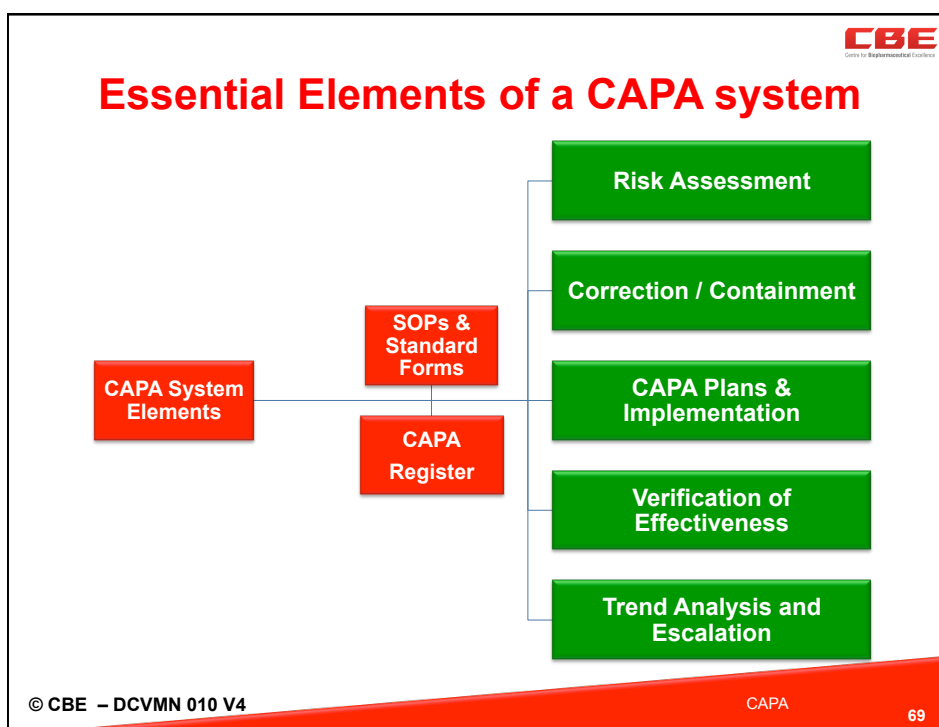

CBE
Center for Biopharmaceutical Excellence

Important “CAPA” Definitions

Correction: Correction refers to repair, rework or adjustment and relates to the disposition of an existing non-conformity, defect, or other undesirable situation



Corrective Action: Action to eliminate the causes of an **existing** non-conformity, defect or other undesirable situation in order to prevent **recurrence**


Preventive Action: Action taken to eliminate the cause of a **potential** non-conformity, defect, or other undesirable situation in order to prevent **occurrence**

Continuous Improvement: Recurring activity to increase the ability to fulfill requirements.


© CBE – DCVMN 010 V4

CAPA

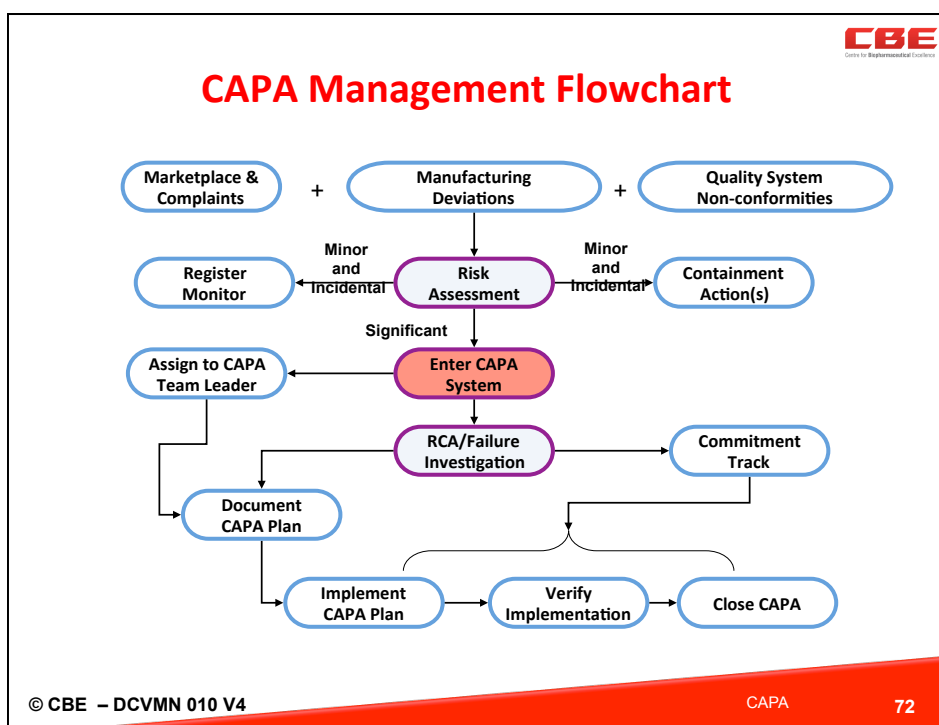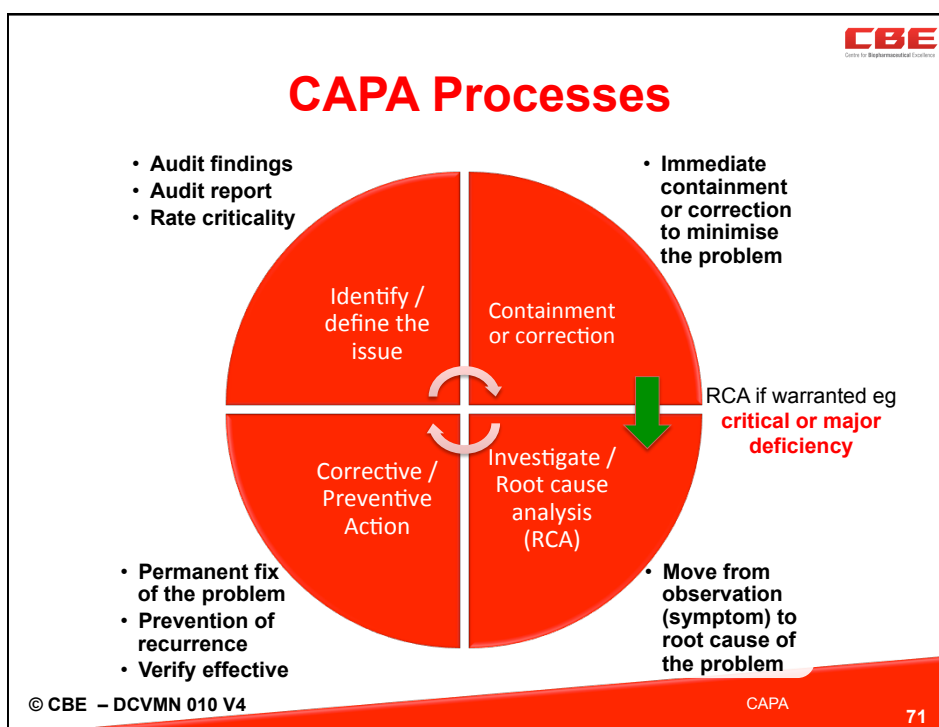




Centre for Biopharmaceutical Excellence



 Centre for Biopharmaceutical Excellence

Steve Williams,
 Director, CBE Pty Ltd
www.cbe-ap.com.au
 ☎ +61(0)417116476
 ✉ steve.williams@cbe-ap.com.au



© CBE – DCVMN 010 V4

73